

Alvi Modified Hohmann Retractor

Designed by Hasham Alvi, MD

Designed for use during minimally invasive anterior hip replacement surgery, the retractor is placed through the capsule, into the femoral head, allowing for retraction of the rectus femoris




The extra bend in the handle allows the assistant to stand on the operative side of the table allowing for ease of handling of the retractor.

PRODUCT NO:

4549

Overall Length: 8.75" (22,2 cm)
Blade Width: .75" (19 mm)

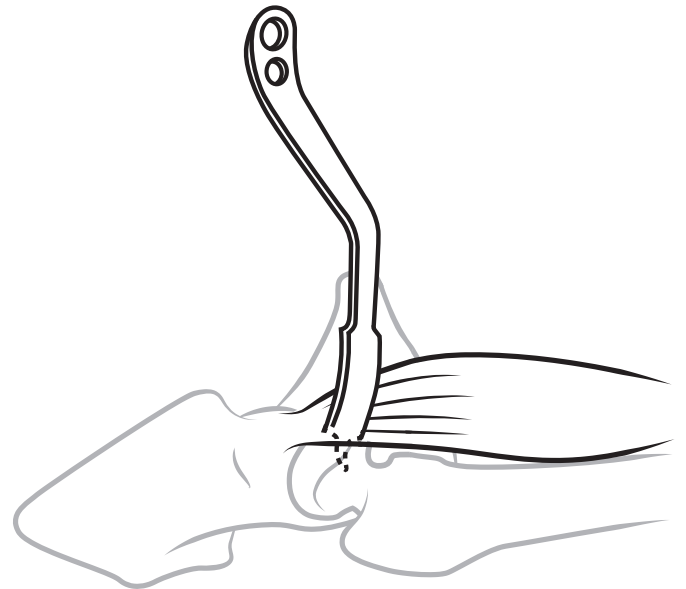

www.innomed.net

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

INNOMED



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362