


Angled Glenoid Retractors – Forked

Designed by R.L. Stowell, MD

Designed to help with exposure of the difficult glenoid and facilitation of glenosphere placement for reverse arthroplasty, with wide and narrow versions to accommodate most glenoid variations

Specifically designed bend allows the retractor to sit along the posterior aspect and deliver the glenoid vault into the surgical field while retracting the humeral head posterior and lateral without excessive force on the glenoid neck, while the working end of the design allows reaming and peripheral bone removal without binding on the reamer.

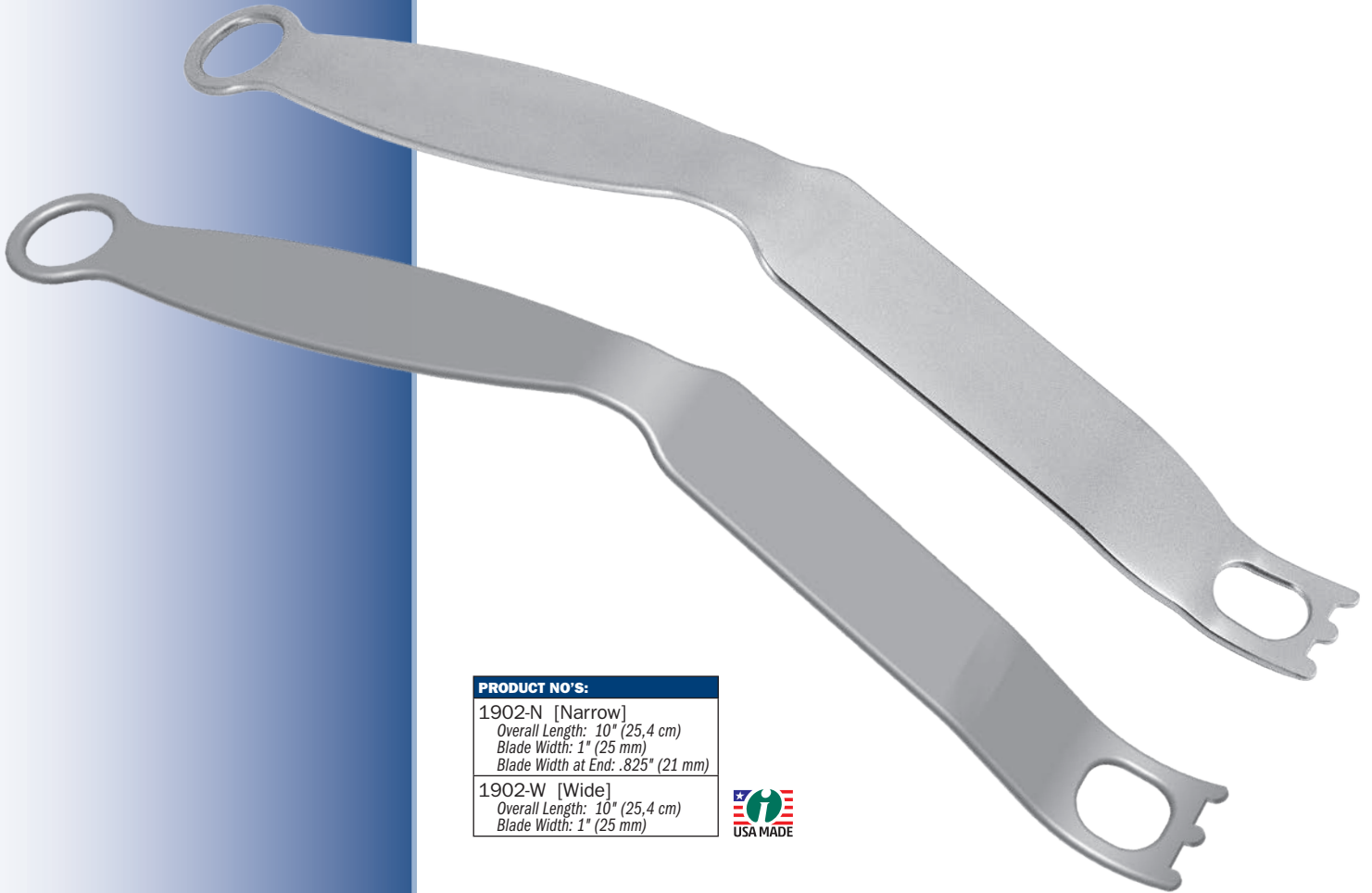

PRODUCT NO'S:
1902-N [Narrow] Overall Length: 10" (25,4 cm) Blade Width: 1" (25 mm) Blade Width at End: .825" (21 mm)
1902-W [Wide] Overall Length: 10" (25,4 cm) Blade Width: 1" (25 mm)



ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

INNOMED



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

1.800.548.2362