


# Bates Needle Holder with Suture Cutter

Designed by James E. Bates, MD



*By trapping the suture and cutting when the forcep is opened, helps to reduce stress on the surgeon's hand*

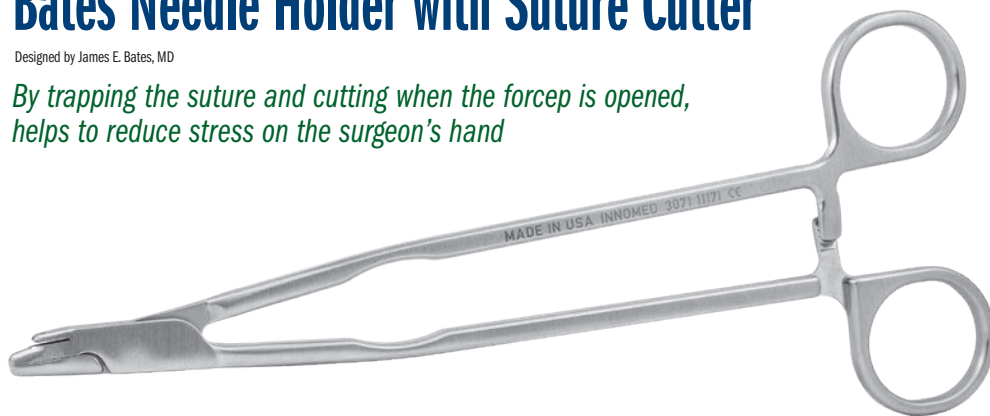

- ▶ No switching between needle driver and scissors, or need for assistant to cut sutures for you
- ▶ Cutting with opening of forceps reduces possibility of damage to surrounding tissues
- ▶ Sliding the instrument down to the suture knot allows quick and consistent 2 mm suture tails
- ▶ Slip the suture strands into the suture cutting slot and slide the closed instrument along until desired length of tail is achieved, then open the instrument to cut the sutures

## PRODUCT NO:

3071  
Overall Length: 8.125" (20,6 cm)  
Jaw Width: .25" (6,4 mm)  
Open Jaw Length: .5" (12,8 mm)

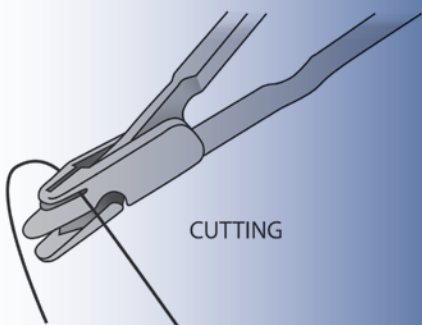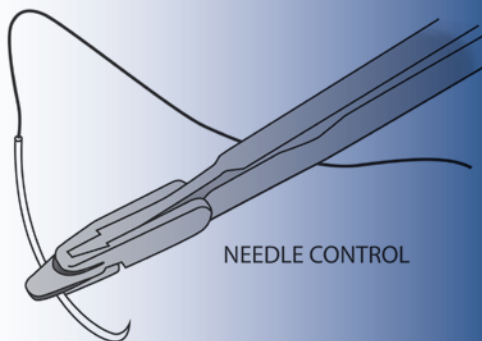

[www.innomed.net](http://www.innomed.net)

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404  
info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74  
info@innomed-europe.com

1.800.548.2362

# INNOMED