

# Cannestra Hip Positioner



Designed by Vincent Cannestra, MD

*Designed for stable positioning of a patient during total hip, MIS and revision surgery*

Low profile allows the surgeon and assistant close proximity to the patient—helpful for MIS surgery.

With the patient laying on the table, the anterior and posterior assemblies can be moved independent of each other. The anterior pad assembly can be adjusted for height. The anterior cushions secure the pelvis against the symphysis pubis and anterior superior iliac spine of the non-operative side. Allows greater flexion of the operative hip.

Attaches directly to the operating table utilizing the included table clamps. The pads are made of semi-dense foam to help prevent pressure points and are sealed with a washable coating. The complete unit is radiolucent and autoclavable except for the foam pads.

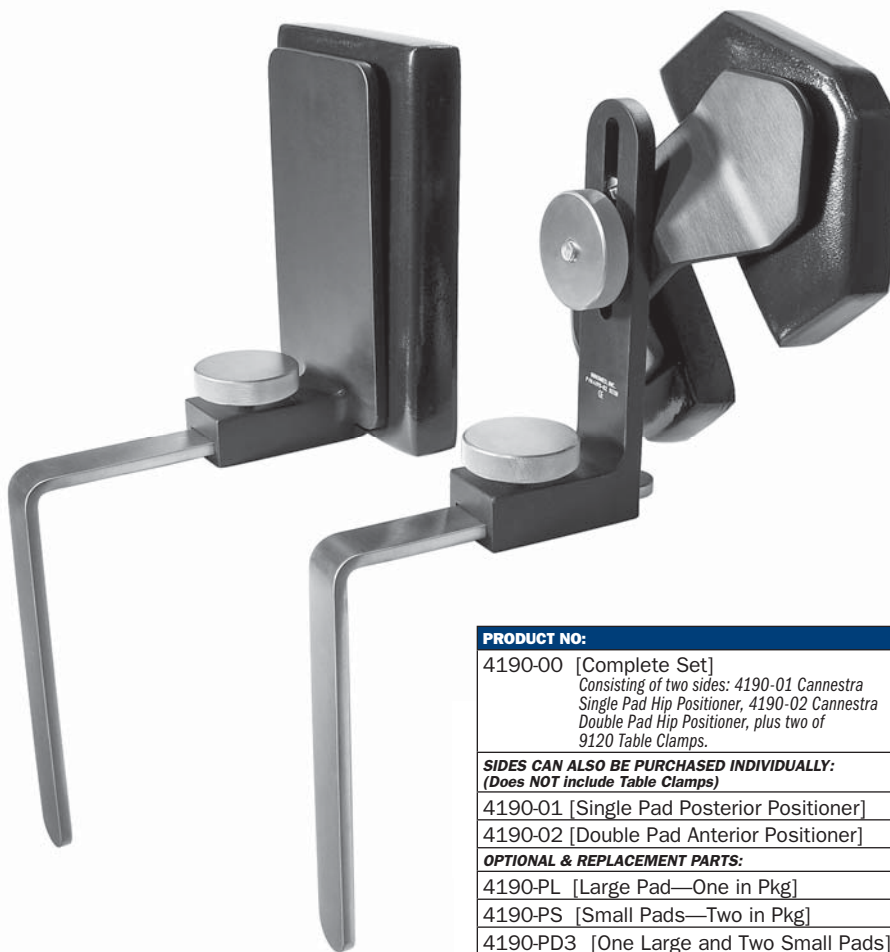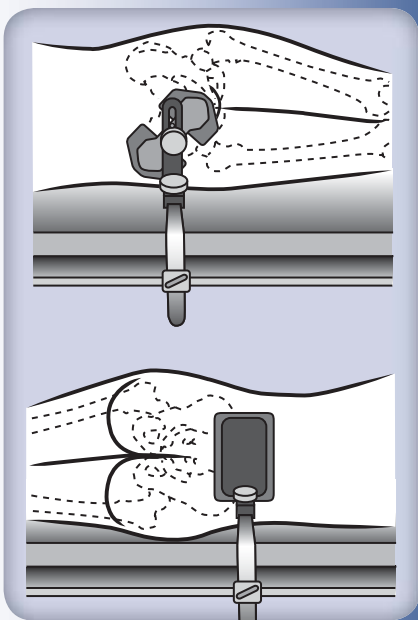

#### PRODUCT NO:

4190-00 [Complete Set]  
Consisting of two sides: 4190-01 Cannestra Single Pad Hip Positioner, 4190-02 Cannestra Double Pad Hip Positioner, plus two of 9120 Table Clamps.

#### SIDES CAN ALSO BE PURCHASED INDIVIDUALLY: (Does NOT include Table Clamps)

4190-01 [Single Pad Posterior Positioner]

4190-02 [Double Pad Anterior Positioner]

#### OPTIONAL & REPLACEMENT PARTS:

4190-PL [Large Pad—One in Pkg]

4190-PS [Small Pads—Two in Pkg]

4190-PD3 [One Large and Two Small Pads]

9120 [Table Clamp only]

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404  
www.innomed.net info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax



# INNOMED

1.800.548.2362