

1-800-548-2362

International (912) 236-0000



www.innomed.net

(912) 236-7766 Fax

©2005 Innomed, Inc., All Rights Reserved.

Innomed

103 Estus Drive, Savannah, GA 31404
E-Mail: info@innomed.net



FREE TRIAL ON MOST INSTRUMENTS
CALL FOR MORE INFORMATION




Femoral Head Removal Pin

Partial threaded pin can be used to help remove a femoral head during total hip surgery. The pin is especially helpful in minimally invasive total hip surgery where access to the femoral head is limited. The pin is attached to a pin driver which clamps onto a Jacob chuck. When the pin is drilled in place, the driver is easily removed from the pin, as the pin is held by a friction ring. The head can be removed by gripping the pin by hand or by using a large pin inserter/extractor. Pin diameter is 5/32".

Product No's:
1310 [Pin]
1205 [Pin Driver]
3030 [Pin Inserter/Extractor]

Used to help remove a femoral head during total hip surgery

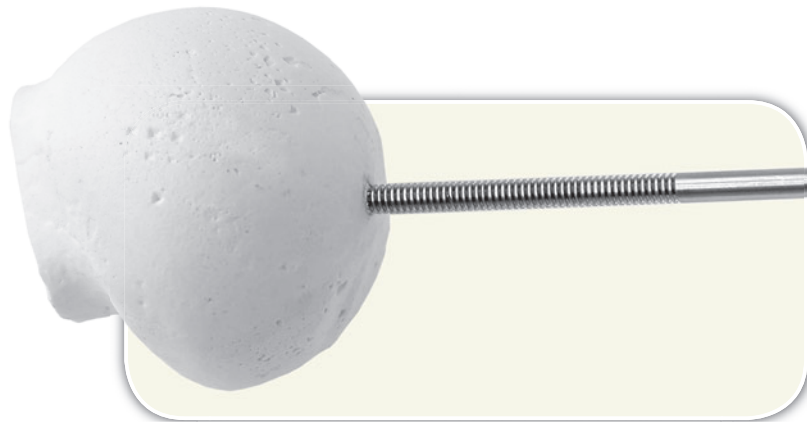

Pin Inserter/Extractor



Pin Driver