

Femoral Head Disengaging Punch

Designed by Brandon Thompson, CST/CFA

Designed to help protect the femoral stem trunnion while removing the femoral head


The delrin pad helps prevent scratching of the femoral stem trunnion.
The punch angle allows for better striking force to help break the taper of the head and stem.



THE DELRIN PAD SHOULD NOT BE USED FOR IMPACTION.

PRODUCT NO:

8626
Overall Length: 9" (22.9 cm)
Shaft Diameter: .5" (12.7 mm)
Punch Platform Offset Angle: 30°
Punch Platform Delrin End: 10 mm x 20 mm

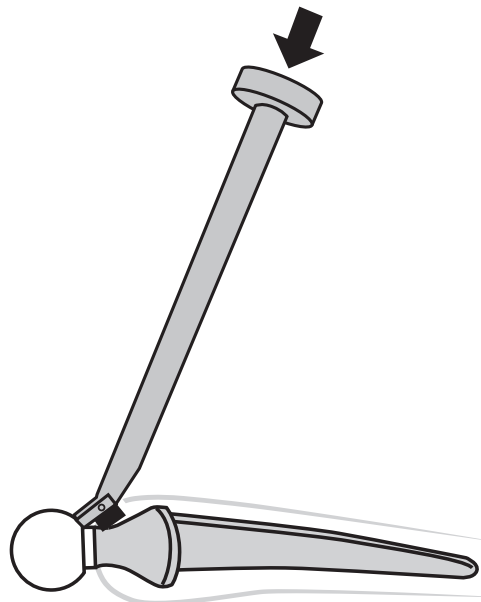

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe

Tel. +41 41 740 67 74
Fax +41 41 740 67 71

1.800.548.2362

INNOMED