

Fromm Femur & Tibia Triangles

Designed by S.E. Fromm, MD *
Extra Small Triangle designed by S.E. Fromm, MD & Kenneth Merriman, MD

Used for femur and tibia positioning during nailing, repairs and fractures

Designed to position and hold the femur and tibia during intramedullary nailing of the tibia, ligament repairs and extremity fractures. Allows knee to be flexed greater than 90° to allow reaming and nail insertion without displacing fracture. The triangles are available in four heights: 8.5", 11", 14", and 16". The three smaller triangles are designed to fit inside the larger triangle for storage. They are supplied with an autoclavable silicone cushioning pad and velcro* straps. The triangles are also radiolucent and gas or steam sterilizable.

*Velcro® is a registered trademark of the Velcro Companies.



PRODUCT NO'S:

2760-00 [Set of 3] Angles: Top 30°, Two Bottom 75°
2760-01 [11"] Base: 6" (15,2 cm), Height: 11" (27,9 cm)
2760-02 [14"] Base: 7" (17,8 cm), Height: 14" (35,6 cm)
2760-03 [16"] Base: 9" (22,9 cm), Height: 16" (40,7 cm)
Sold Separately - Not In Set:
2760-XS [8.5"] Base 5" (12,7 cm), Height: 8.5" (21,6 cm)

Replacement Parts:

2760-P [Silicone Pad]
2760-S [Straps] Package of 18 - 6 Blue / 12 Green
8100-P [Green Straps for Femur, Long] Pkg of 10
8120-P [Blue Straps for Tibia, Short] Package of 10
8120-SP [Straps for 2760-XS] Package of 10



— Over —

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

INNOMED



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net


912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71


1.800.548.2362




Fromm Femur & Tibia Triangles

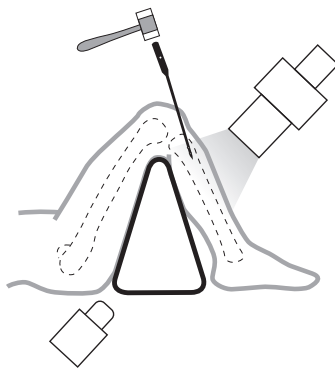


Tibial Nailing

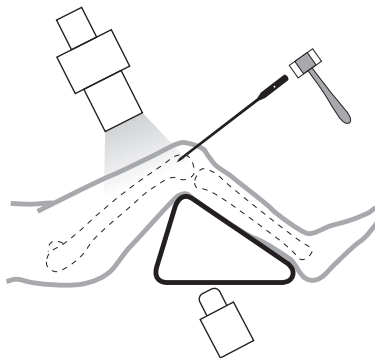

Retrograde Femoral Nailing

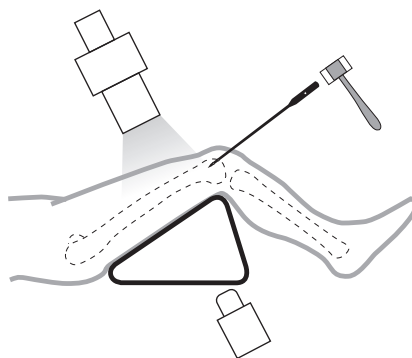

Retrograde Femoral Nailing
Triangle holds femur reduced
(prevents sagging)

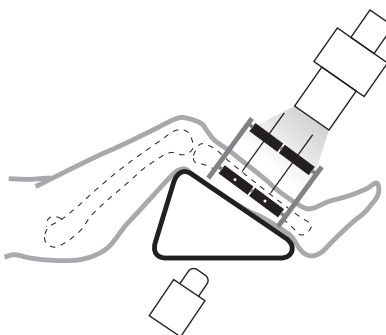

Triangle Holds Tibia Reduced For:

- Open Reduction and Internal Fixation (ORIF) of tibia
- Application of uniplane or multiplane external fixator
- Knee ligament repairs and/or reconstruction

— Over —

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe

Tel. +41 41 740 67 74
Fax +41 41 740 67 71

1.800.548.2362

INNOMED