


Lachiewicz-Hoover Patella Retractor

Helps to retract and prevent harm to the polyethylene implant from scratches, gouges, and bovie cautery during debridement and removal of other components

In primary TKA, it is used after implantation of components so the surgeon can debride the lateral gutter of synovium and for exposure for hemostasis laterally. Also extremely valuable for revisions in which a well fixed patella component is to be retained.

Designed By
Paul Lachiewicz, MD &
Steve Hoover, MD

Product No: 3235
Overall Length: 6.5"
Patella Pad Diameter: 1.5"
Patella Pad Thickness: .25"

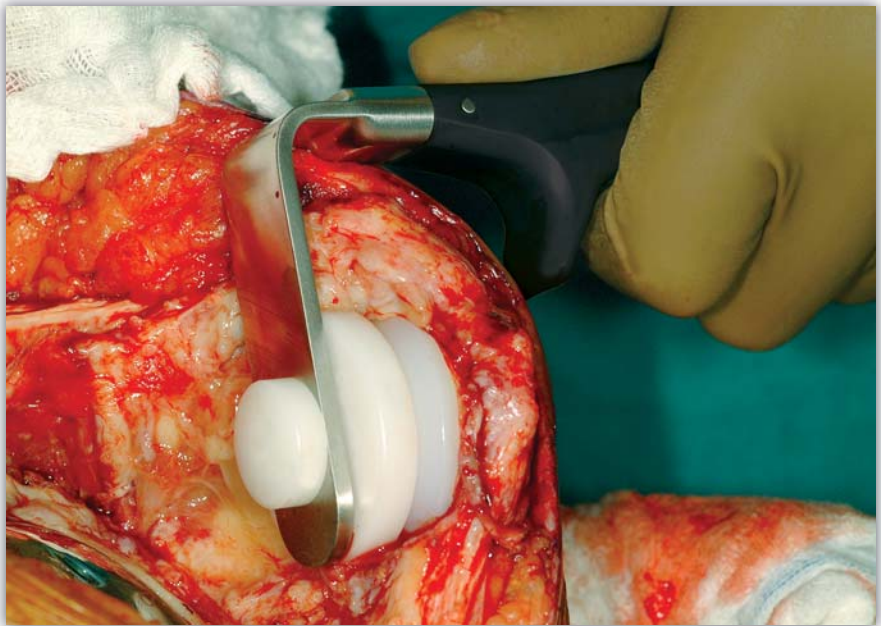

FREE TRIAL
ON MOST INSTRUMENTS
CALL FOR MORE INFORMATION



©2008 Innomed, Inc., All Rights Reserved

Phone **800 548 2362**
912 236 0000
Fax 912 236 7766
Internet www.innomed.net
eMail info@innomed.net

Innomed
ORTHOPEDIC INSTRUMENTS
103 Estus Drive Savannah GA 31404

