

Nicholson Small Bone and Shoulder Cement Removal Instruments

Designed by Gregory Nicholson, MD

Designed to facilitate cement removal in smaller diameter bone of the humerus, ulna, and smaller implant geometries




PRODUCT NO'S:

| |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Gouges Overall Length: 9" (22,9 cm) Gouges Handle Length: 4" (10,2 cm) |
| 5251-00 [Complete Set with Case] |
| Also Available Individually: |
| 5251-05 [Extra Small Gouge] Gouge Width: 5 mm Chisel Depth: .675" (1,6 cm) |
| 5251-07 [Small Gouge] Gouge Width: 7 mm Chisel Depth: 1" (2,5 cm) |
| 5251-09 [Medium Gouge] Gouge Width: 9 mm Chisel Depth: 1" (2,5 cm) |
| 5251-11 [Large Gouge] Gouge Width: 11 mm |
| 5252-07 [Small Gouge w/Splitter] Gouge Width: 7 mm Splitter Height: 4 mm Chisel Depth: 1" (2,5 cm) |
| 5252-09 [Medium Gouge w/Splitter] Gouge Width: 9 mm Splitter Height: 5 mm Chisel Depth: 1" (2,5 cm) |
| 5252-11 [Large Gouge w/Splitter] Gouge Width: 11 mm Splitter Height: 6 mm Chisel Depth: 1" (2,5 cm) |
| 5254 [Backhook] Overall Length: 12.5" (31,8 cm) Handle Length: 4.5" (11,4 cm) Shaft Diameter: 4 mm |
| 5255 [Footed Impactor] Foot Pad Size: 8.5 mm x 11.5 mm Shaft Diameter: 8.5 mm (21,6 cm) Overall Length: 12.75" (32,4 cm) Handle Length: 4.5" (11,4 cm) |
| 5253 [Case for Set] |

- ▶ Reverse bevel tip helps the gouge to slide between the bone and cement.
- ▶ T-shaped Gouge-Splitter allows the gouge to slide between the cement and bone and vertically split the cement mantle to facilitate removal.
- ▶ Small diameter widths and curvatures more closely match shoulder and elbow implants and smaller bone diameters.
- ▶ Shorter length allows for better control and access.



Used to help remove a humeral implant by impacting the medial collar of the prosthesis – helps provide a very direct parallel force to the implant for removal

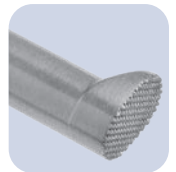

www.innomed.net

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

INNOMED



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362

© 2023 Innomed, Inc.