
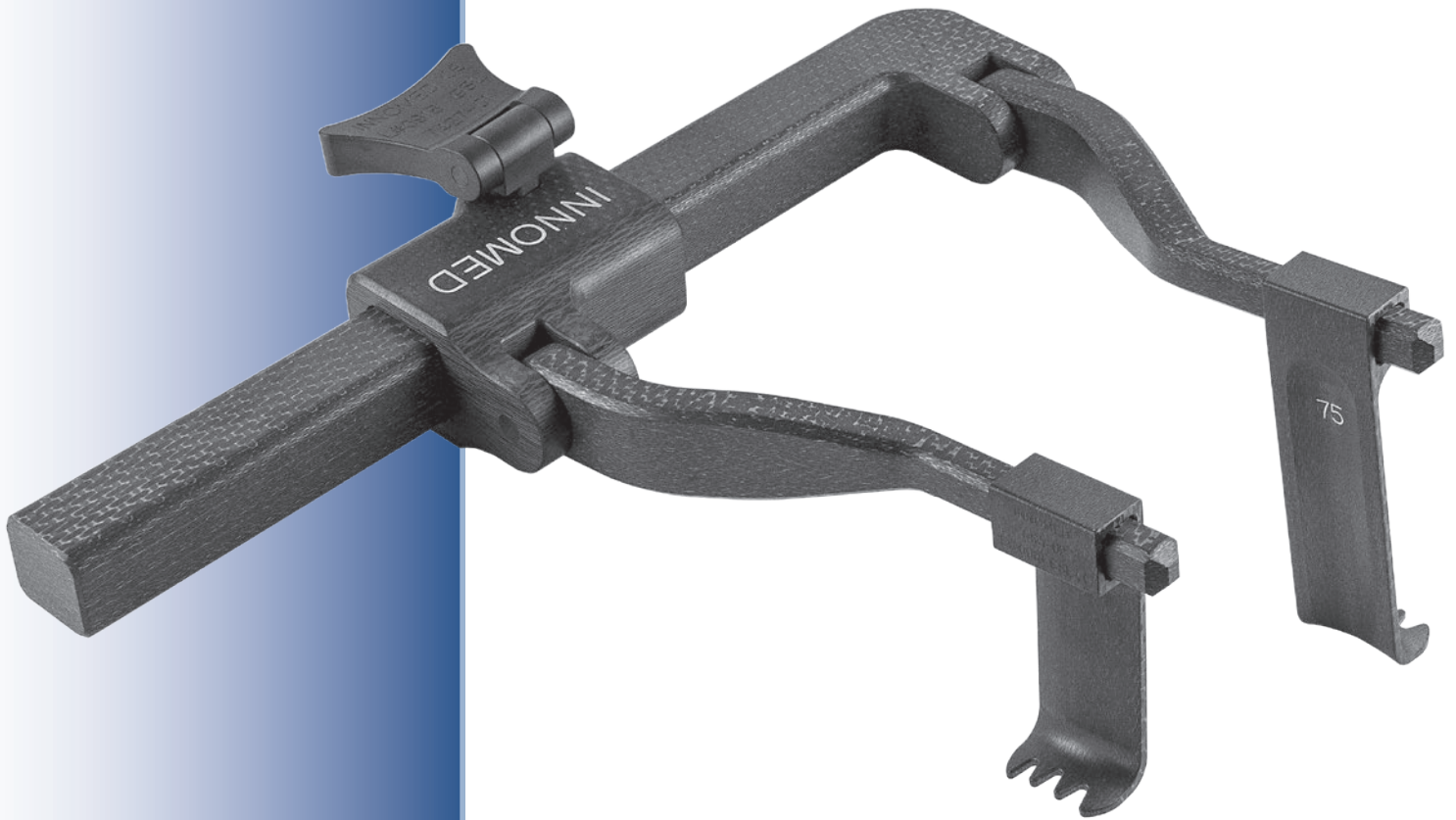


OrthoLucent™ Self-Retaining Ratchet Frame Assembly

Designed for wound exposure, the 100% radiolucent retractor frame, arms, and blades can be kept in place while using image intensification or taking an x-ray

Made of strong, lightweight carbon fiber PEEK composite material, which is completely radiolucent, helps to prevent from marring component surfaces, and can be steam sterilized.

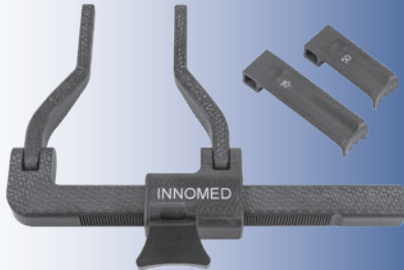

- ▶ Two blade depths: 50 & 75 mm
- ▶ Arms rotate 180°
- ▶ Ratchet Knob pulls out slightly for turning, and pushes in for locking
- ▶ Knob and Arm Unit are removable for cleaning

Ratchet Knob **Pulled Out:**
Turnable and Unlocked



Ratchet Knob **Pushed In:**
Locked and Unturnable



PRODUCT NO'S:	
7427-01 [Frame] (Includes Handle)	Dimensions (Flat): 9.625" x 7.25" (24,5 cm x 18,4 cm)
7427-02 [50 mm Blade]	Dimensions: 50 mm Deep X 25 mm Wide
7427-03 [75 mm Blade]	Dimensions: 75 mm Deep X 25 mm Wide

MADE EXCLUSIVELY
FOR INNOMED IN
SWITZERLAND

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

Scan to
Launch Our
Website



INNOMED



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

1.800.548.2362