

Percutaneous Achilles Repair Forceps

FOR LIMITED OPEN ACHILLES TENDON REPAIR

Designed by James A. Amis, MD

Designed to help improve accuracy during percutaneous repair of Achilles tendon ruptures



The bump on the lateral side of each loop allows the surgeon to palpate the exact center of the loop, proximal to distal, and drop a needle just below (patient is prone) or anterior to the bump for the starting point, and aim to just below the bump on the opposite side.



Lateral Bump

The bump on the lateral side of each loop allows the surgeon to palpate the exact center of the loop, proximal to distal, and drop a needle just below (patient is prone) or anterior to the bump for the starting point, and aim to just below the bump on the opposite side

PRODUCT NO:
8235
Overall Length: 9.625" (24,4 cm)

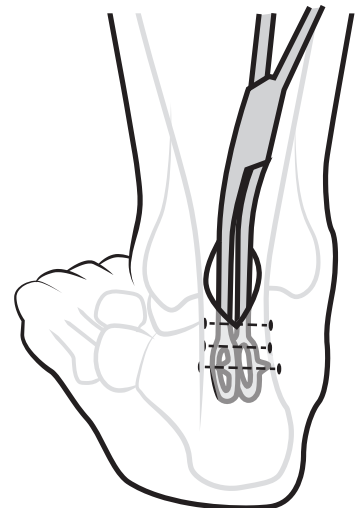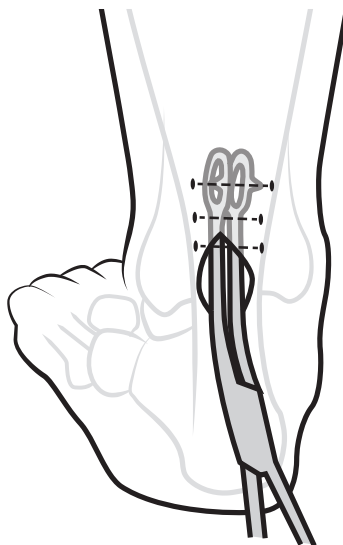

ISO 9001:2008 • ISO 13485:2003

FREE TRIAL ON MOST INSTRUMENTS

Scan to
Launch Our
Website



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71



1.800.548.2362

© 2016 Innomed, Inc.