


Wilson Condylar Gauge

Designed by Ralph Wilson, MD

Designed to measure the posterior femoral condyle after the posterior cuts have been made in total knee arthroplasty

By measuring the depth of the residual condyle, the surgeon can resect excessive bone and measure the bone remaining to avoid impingement of the condyle against the tibial component which could impair knee flexion. The gauge is applied to the inferior or posterior cut surface of the femoral condyle, and the back to front residual bone is measured and then removed as needed. Measures to 30 mm depth.

PRODUCT NO:
1194
Overall Length: 6" (15,2 cm)
Width: .568" (14,4 mm)

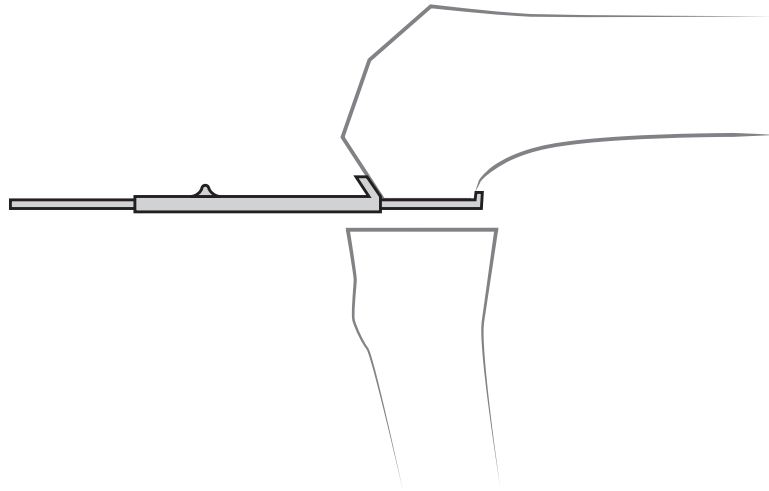

ISO 9001:2008 • ISO 13485:2003

FREE TRIAL ON MOST INSTRUMENTS

Scan to
Launch Our
Website



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

INNOMED

1.800.548.2362