
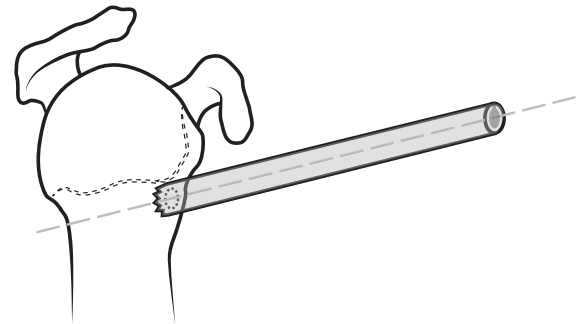


Argintar Bicep Tenodesis Sleeve

Designed by Evan Argintar, MD

Designed to help facilitate mini-open sub-pectoral bicep tenodesis—by maintaining the trajectory of the drill with the serrated end of the sleeve, the drilled humeral holes are easily found with standard percutaneous placement of the bicortical button

Once flipped, the slotted cut out in the sleeves makes detachment of the button applicator possible, helping with efficient and reproducible mini-open bicep tenodesis using button technique.

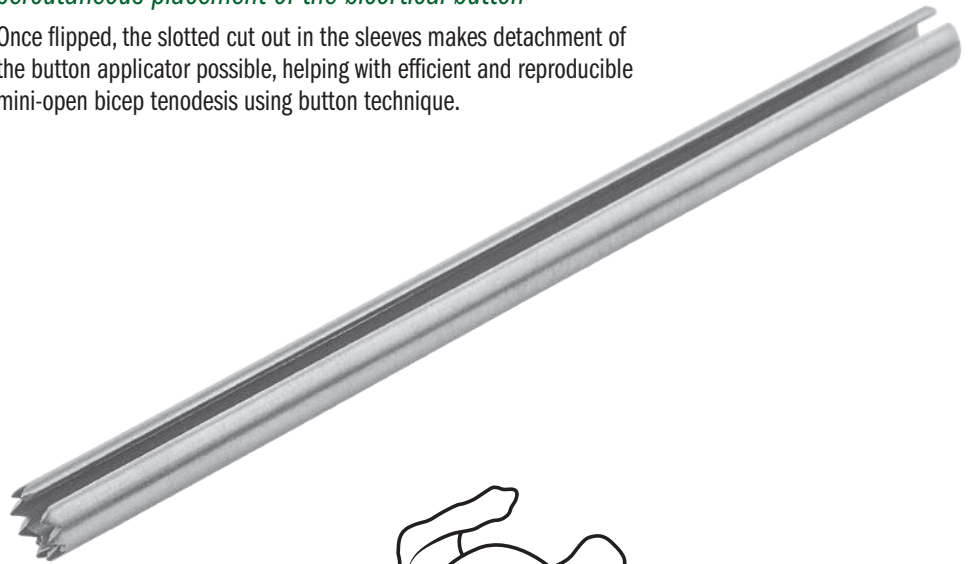

PRODUCT NO:

5835
Overall Length: 4" (10,2 cm)
External Diameter: 6,35 mm
Internal Diameter: 5 mm
Slot: 2,75 mm



ISO 9001:2008 • ISO 13485:2003

FREE TRIAL ON MOST INSTRUMENTS

Scan to
Launch Our
Website



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

INNOMED

1.800.548.2362