




Belfast Sagittal Plane Positioner

Designed by David Beverland, FRCS


A sturdy and stable patient support system for posterior approach total hip arthroplasty in the lateral decubitus position

- ▶ Very secure and easy to tighten
- ▶ Accommodates the very obese patient
- ▶ Does not attach to the table, making it compatible with all OR tables

- ▶ Anterior ASIS support unit allows both height and depth adjustment, eliminating the need to adjust the baseplate locking mechanism once secured. The depth adjustment allows extension of up to 150 mm beyond the support base, providing space to avoid the abdominal apron.
- ▶ Posterior support unit offers both Sacral and Iliac Crest support, connected together as one unit, and provides both height and rotation adjustments, as well as a hinged pivot where the iliac crest side attaches the unit to the post, that can be loosened and adjusted to allow rotation in 5° increments.
- ▶ Allows access to the PSIS before and during surgery. This can help to ensure that the sagittal plane of the pelvis is horizontal at the time of cup implantation.

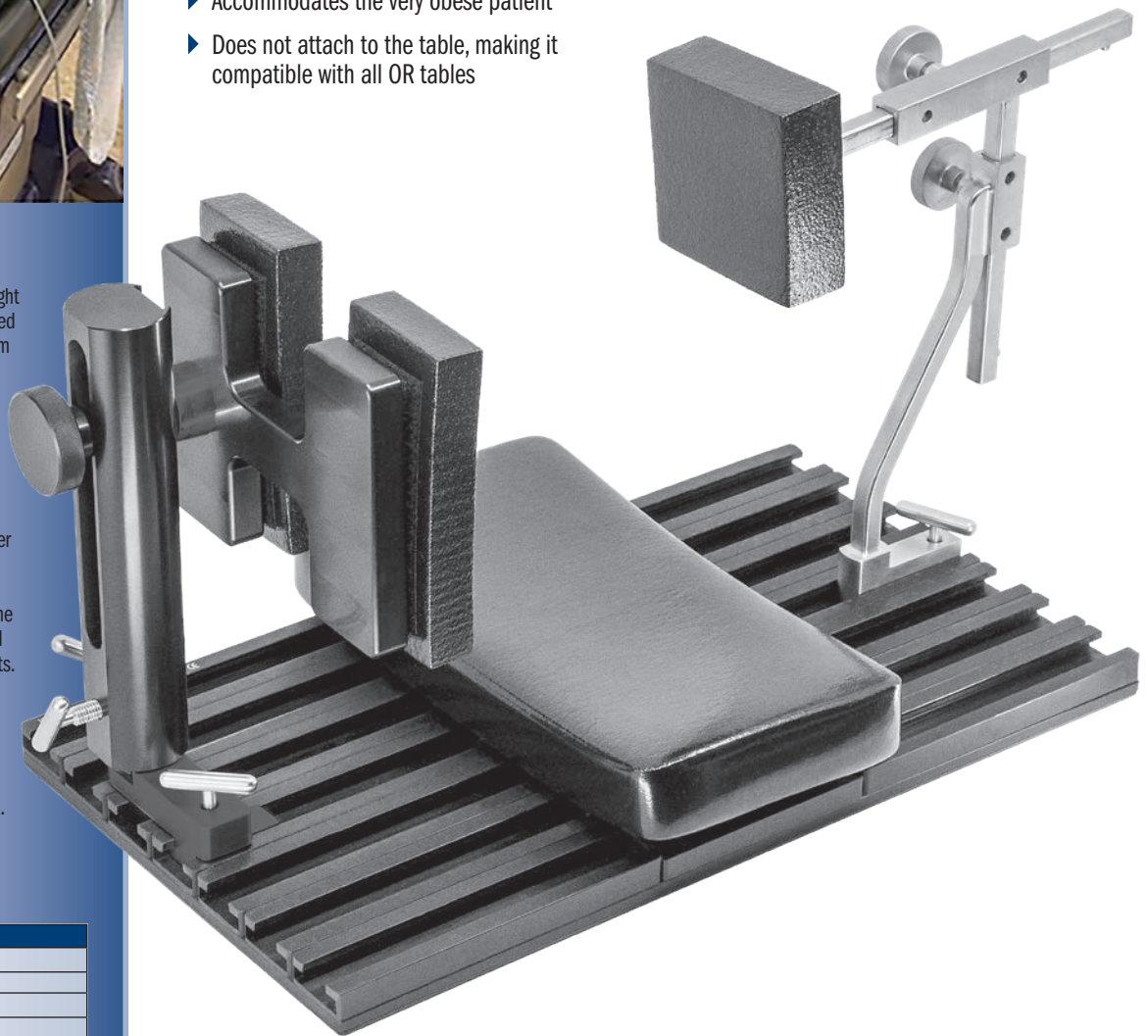

PRODUCT NO'S:	
4170-00	[Complete Set]
Set Includes / Available Individually:	
4170-03	[Anterior Upright Support]
4170-04	[Anterior Plane Support]
4170-05	[Anterior Clamp Support]
4170-AKS	[Anterior Knob Screw]
<i>Two (2) included in Set, One (1) with this product number</i>	
4170-AP	[Anterior Plane Pad]
4170-06	[Posterior Sagittal Plane Support]
4170-07	[Posterior 9.5" Post]
4170-PKS	[Posterior Knob Screw]
4170-08	[Posterior Base]
4170-09	[Posterior Angle Adjuster]
4170-T	[Posterior T-Handle Screw]
4170-PP	[Posterior Support Pad]
<i>Two (2) included in Set, One (1) with this product number</i>	
4150-PS	[Post Screw]
<i>Three (3) included in Set, One (1) with this product number</i>	
4050-BP	[20" Baseplate Only]
4050-LPD	[Hip Positioner Large Pad]

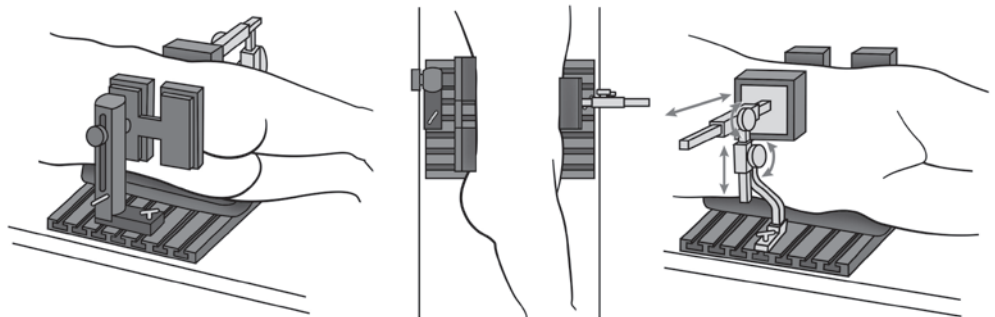

www.innomed.net

FREE TRIAL ON MOST INSTRUMENTS

INNOMED

ISO 13485:2016



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362

© 2023 Innomed, Inc.