


Bicos Meniscal Repair Retractor

Designed by James Bicos, MD

A popliteal retractor specifically designed for meniscal repair or access to the posterior knee

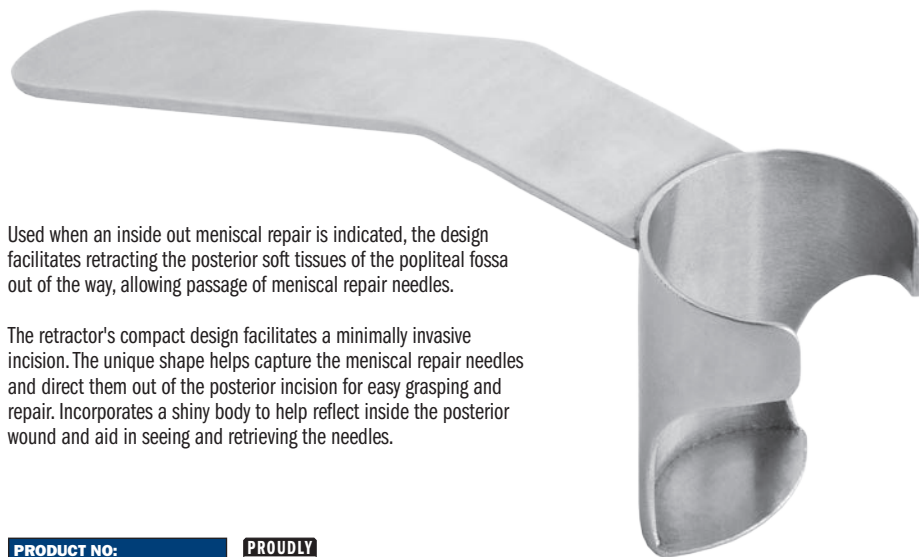

Used when an inside out meniscal repair is indicated, the design facilitates retracting the posterior soft tissues of the popliteal fossa out of the way, allowing passage of meniscal repair needles.

The retractor's compact design facilitates a minimally invasive incision. The unique shape helps capture the meniscal repair needles and direct them out of the posterior incision for easy grasping and repair. Incorporates a shiny body to help reflect inside the posterior wound and aid in seeing and retrieving the needles.

PRODUCT NO:

2731

Overall Length: 5" (12,7 cm)

Depth: 1.625" (4,1 cm)

Diameter: 28 mm

**PROUDLY
MADE
IN THE
USA**

ISO 9001:2008 • ISO 13485:2003

FREE TRIAL ON MOST INSTRUMENTS

Scan to
Launch Our
Website



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

INNOMED

1.800.548.2362