


Desai Curette Osteotomes


Designed by Sarang Desai, DO

Designed to remove bone and cartilage, helpful for preparing joint surfaces for fusion, allowing easy removal of osteophytes and cartilage without having to switch instruments

The osteotome portion also can be used to "feather" the subchondral surface to expose bleeding bone. It is also useful in instances of obtaining autograft, as it can be used to create a bone window and then remove cancellous bone.

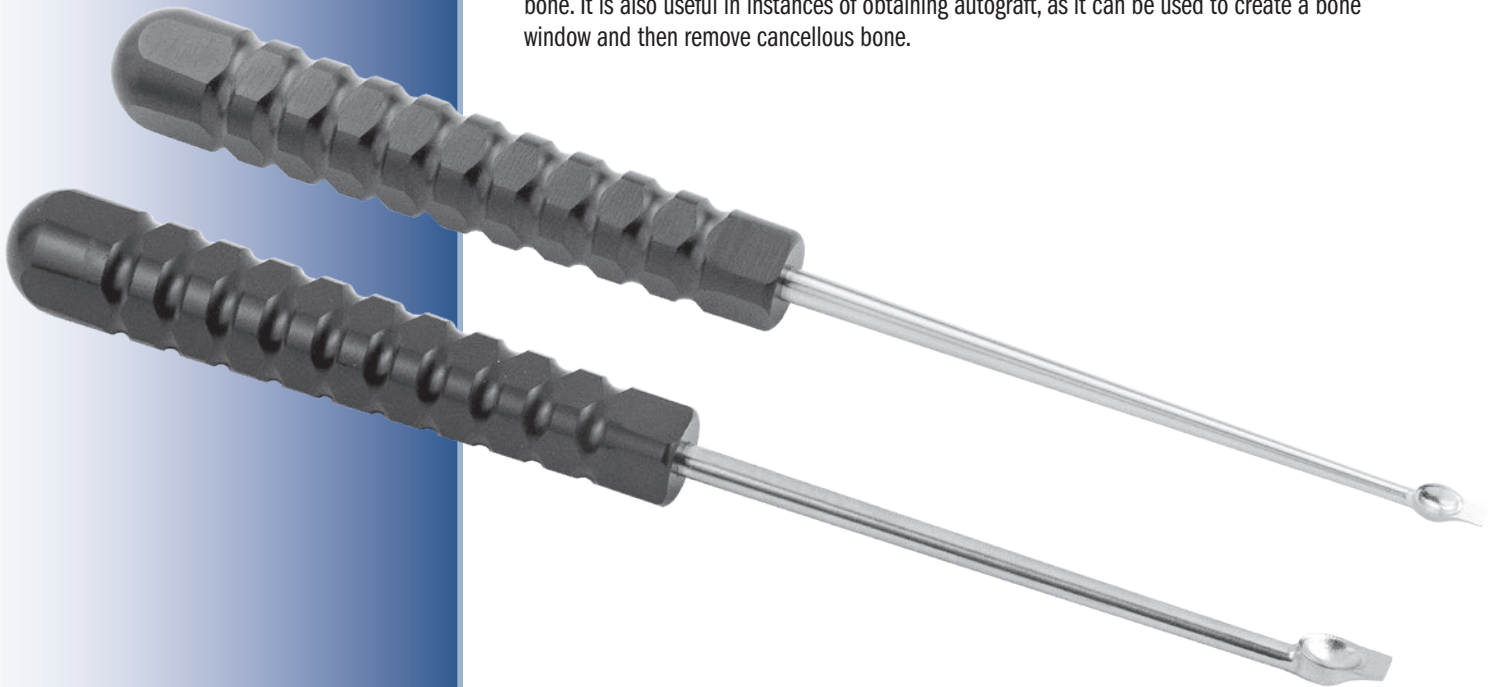

PRODUCT NO'S:	
5241	[5 mm]
Overall Length: 8.25" (21 cm)	
Handle Length: 4.25" (10,8 mm)	
Cup: 5 x 6 mm	
Osteotome Width: 3.5 mm	
Osteotome Length: 3.5 mm from edge of cup	
5242	[8 mm]
Overall Length: 8.25" (21 cm)	
Handle Length: 4.25" (10,8 mm)	
Cup: 8 x 10 mm	
Osteotome Width: 6.5 mm	
Osteotome Length: 3 mm from edge of cup	

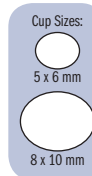

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

Scan to
Launch Our
Website



INNOMED 



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

1.800.548.2362