

## Ludloff/Mau Osteotomy Fixation Clamp

Designed by A. Austin

*Used after lateral hallux valgus correction of the metatarsal, the clamp allows for osteotomy fixation and cannulated screw guide wire direction*

Clamp fixates the osteotomy to hold the correction, and the 15° slanted cannulated k-wire guide allows the surgeon to place the guide wire for the cannulated screw perpendicular to the osteotomy for final fixation of the osteotomy.


**PRODUCT NO:**

1812

Cannula Accepts K-wire up to: .045" (1,1 mm)

Overall Length: 5" (12,7 cm)

MADE EXCLUSIVELY  
FOR INNOMED IN  
GERMANY

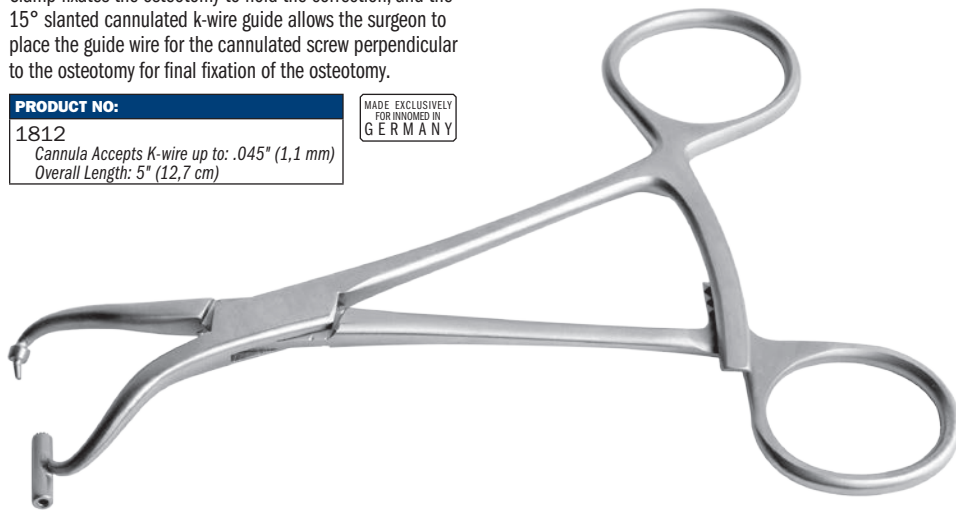

ISO 9001:2008 • ISO 13485:2003

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404  
www.innomed.net info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax

Scan to  
Launch Our  
Website



Innomed-Europe Tel. +41 41 740 67 74  
Fax +41 41 740 67 71

**INNOMED**

1.800.548.2362