


Mantis Screwdriver Distractor

Designed by J. Albert Diaz, MD

*Designed to help provide stable distraction across difficult-to-reduce fractures using two seated screwdrivers**

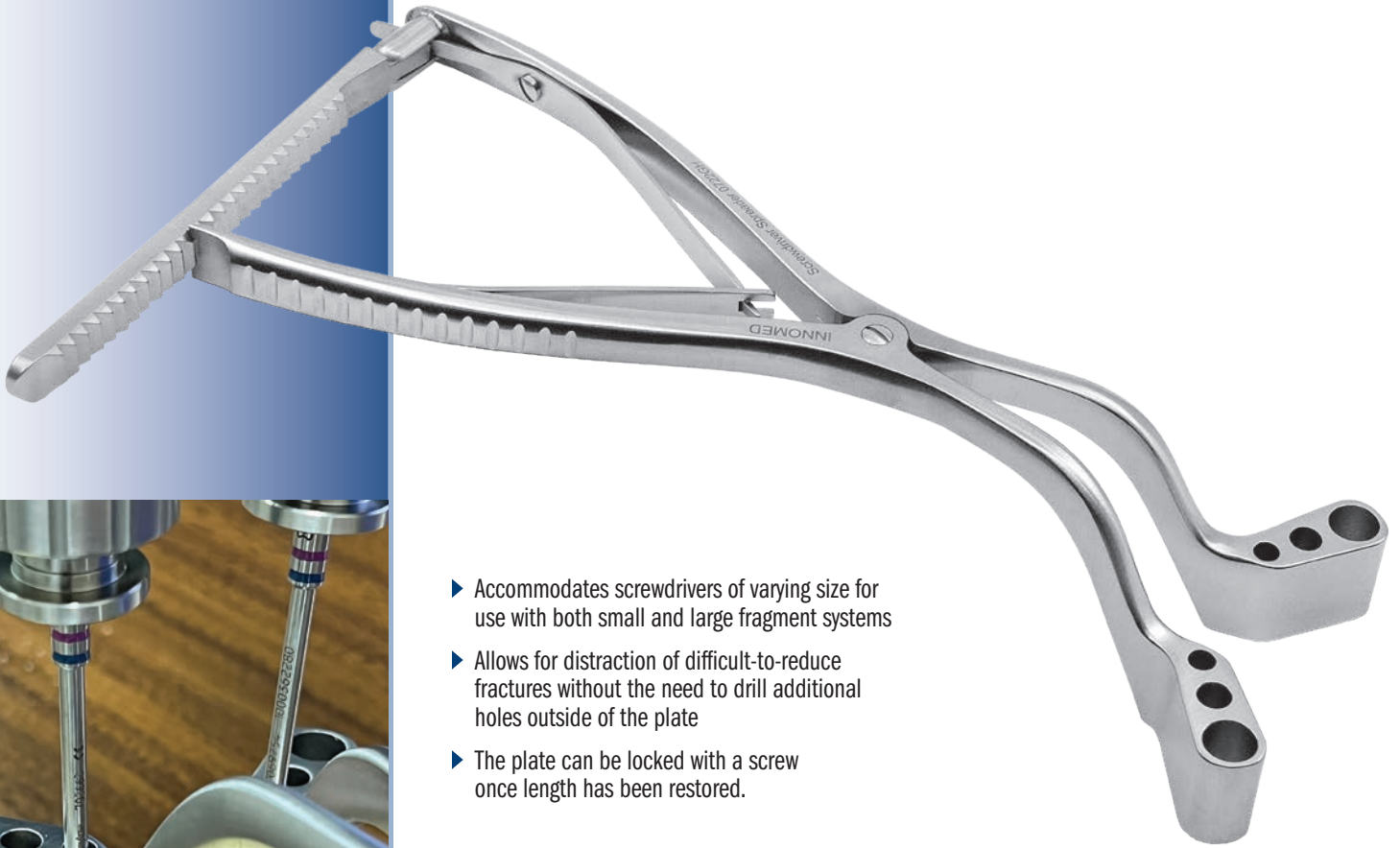

- ▶ Accommodates screwdrivers of varying size for use with both small and large fragment systems
- ▶ Allows for distraction of difficult-to-reduce fractures without the need to drill additional holes outside of the plate
- ▶ The plate can be locked with a screw once length has been restored.

*Screwdrivers not included.

PRODUCT NO:

3654
Overall Length: 7.5" (19,1 cm)
Pin Hole Diameters: 4.5, 5.5, & 8.5 mm
Leg & Pin Hole Depth: .7" (17,5 mm)

MADE EXCLUSIVELY
FOR INNOMED IN
GERMANY



ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS



103 Estus Drive, Savannah, GA 31404
www.innomed.net info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
Fax +41 41 740 67 71

1.800.548.2362

INNOMED 