

# Meckel Posterior Stabilizing Knee Retractor

Designed by Christopher M. Meckel, MD

*Designed to provide enhanced anterior translation of the tibia when doing posterior stabilized total knee replacement*

The 15 mm deep blade section of the retractor is used to lever the tibia forward (by resting the tip on the posterior tibia and the middle blade section block levering off the distal femur) after the box cut has been made in the distal femur.


**PRODUCT NO:**

4538

Blade Width: 20 mm

Blade Depth: 15 mm

Overall Length: 10" (24,5 cm)

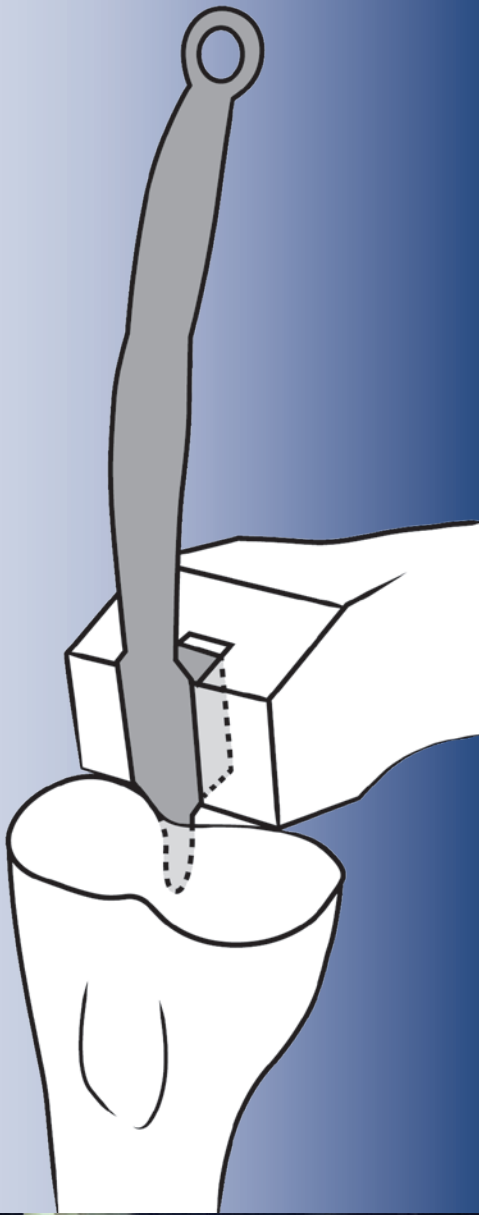

ISO 13485:2016

[www.innomed.net](http://www.innomed.net) FREE TRIAL ON MOST INSTRUMENTS

# INNOMED



103 Estus Drive, Savannah, GA 31404  
info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74  
info@innomed-europe.com

1.800.548.2362