


Multi-Purpose Hip & Knee Retractors

Designed by Vasilius Mathews, MD

Designed for use in both hip and knee arthroplasty procedures

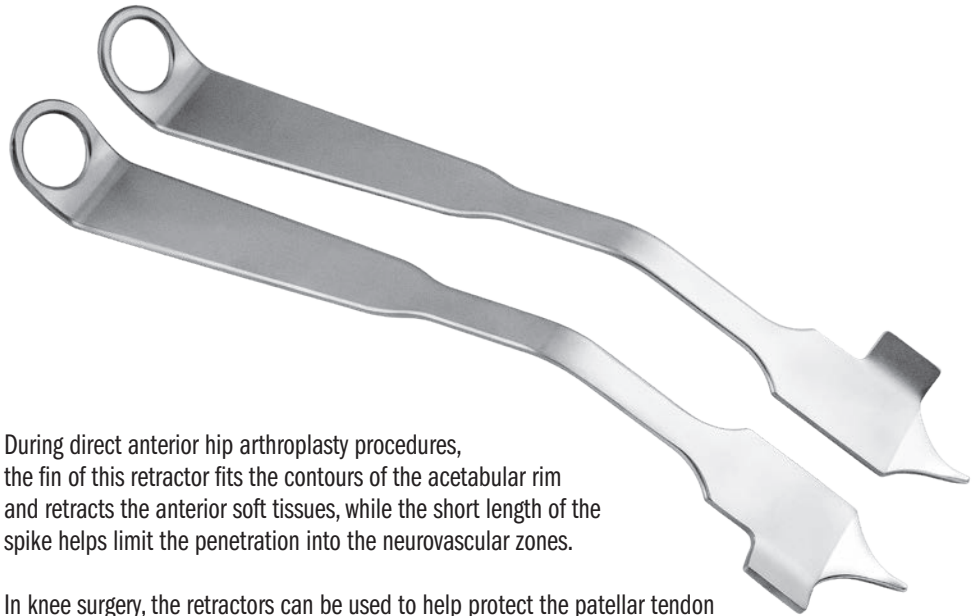


During direct anterior hip arthroplasty procedures, the fin of this retractor fits the contours of the acetabular rim and retracts the anterior soft tissues, while the short length of the spike helps limit the penetration into the neurovascular zones.

In knee surgery, the retractors can be used to help protect the patellar tendon behind the fin at the lateral tibial border. Also useful as a soft-tissue and fat pad retractor during prosthesis implantation, helping to ensure a dry cancellous bed for cementation.

PRODUCT NO'S:

4554-L [Left]
Overall Length: 11.25" (28,6 cm)
Blade Width: 1.5" (38 mm)

4554-R [Right]
Overall Length: 11.25" (28,6 cm)
Blade Width: 1.5" (38 mm)

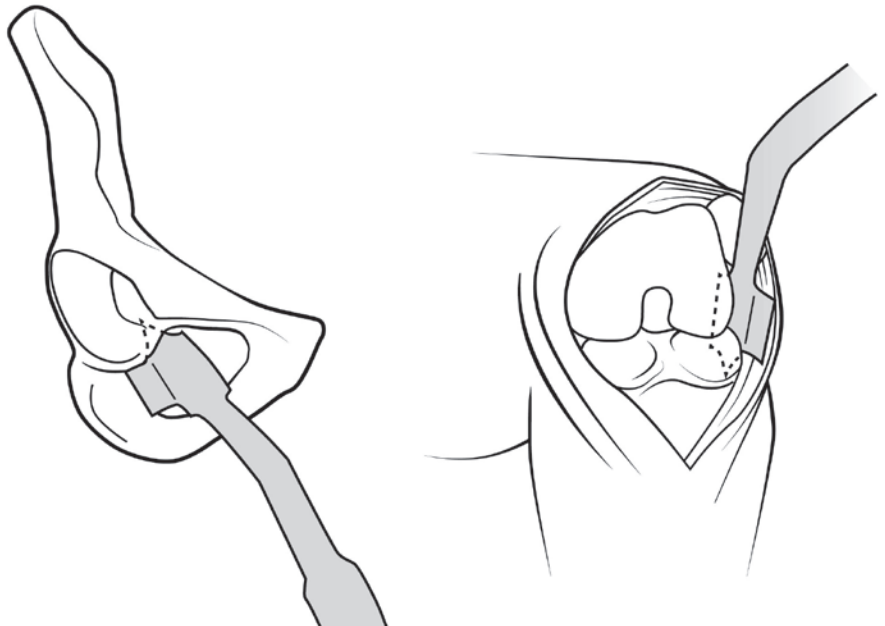

www.innomed.net



FREE TRIAL ON MOST INSTRUMENTS

INNOMED

ISO 13485:2016



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362

© 2024 Innomed, Inc.