
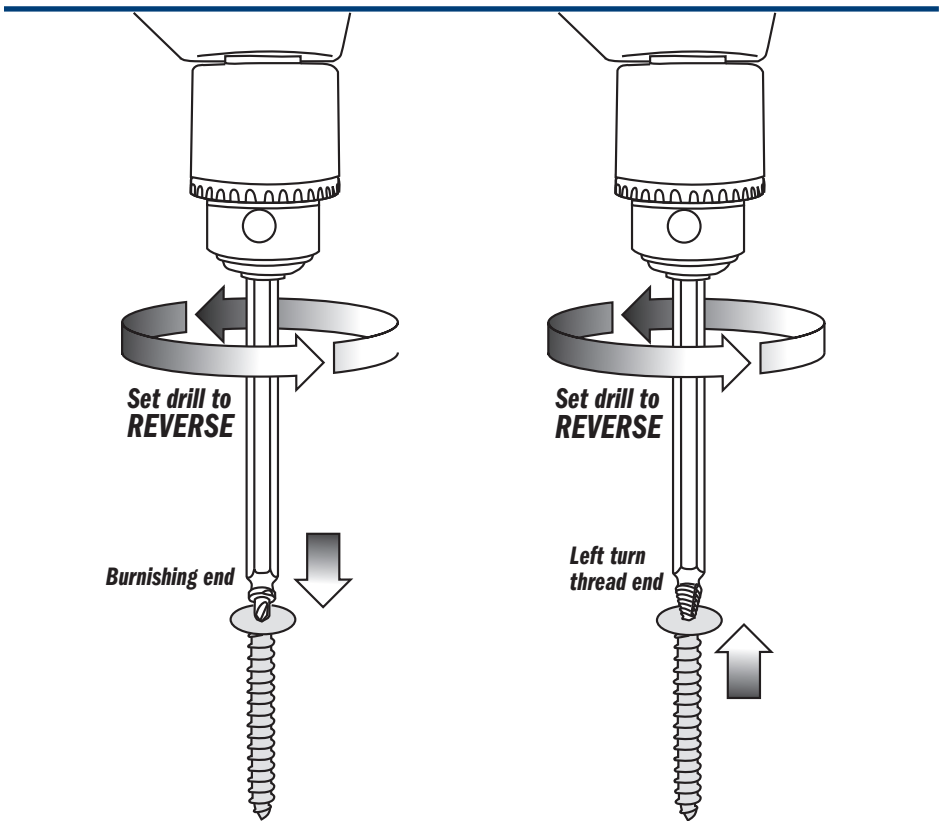


# Screw Extractors

*Designed to help remove screws with stripped or damaged heads*

- ▶ Extractors must be used with drill in reverse.
- ▶ Screw head is reamed with burnishing end, and is then removed with the left turn thread end.
- ▶ Care must be taken to burnish no more than 1/16" (1.5 mm) deep, as burnishing too deep can weaken the screw head.

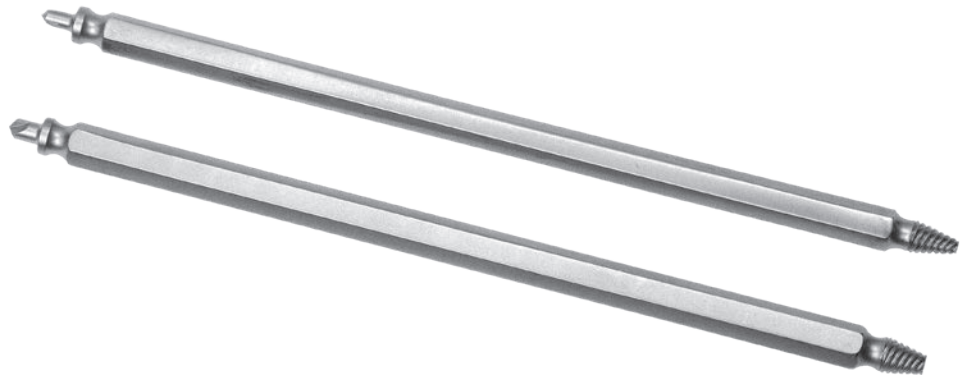

Burnishing Ends

Left Turn Thread Ends

**PRODUCT NO'S:**

7250-01 [2.5 mm Screw Extractor]  
Overall Length: 6" (15,2 cm)

7250-02 [3.5 mm Screw Extractor]  
Overall Length: 6" (15,2 cm)



ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

**INNOMED** 



103 Estus Drive, Savannah, GA 31404  
www.innomed.net info@innomed.net

912.236.0000 Phone  
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74  
Fax +41 41 740 67 71

1.800.548.2362