

Vosburg Cannulated Periarticular Clamp

Designed by Caleb Vosburg, MD

Cannulated clamp tips allow passage of pins up to 7/64" (2.8 mm)

By compressing the fracture with the clamp and then passing two pins, the clamp can then be removed to allow more working room and versatility when applying a plate.

PRODUCT NO:

1864


Overall Length: 13" (33 cm)

Handle Length: 8" (20,3 cm)

Ratcheted Opening from 2" to 3.5" (5,4 to 7,6 cm)

Accepts Pins up to: 7/64" (2.8 mm)

MADE EXCLUSIVELY
FOR INNOMED IN
GERMANY

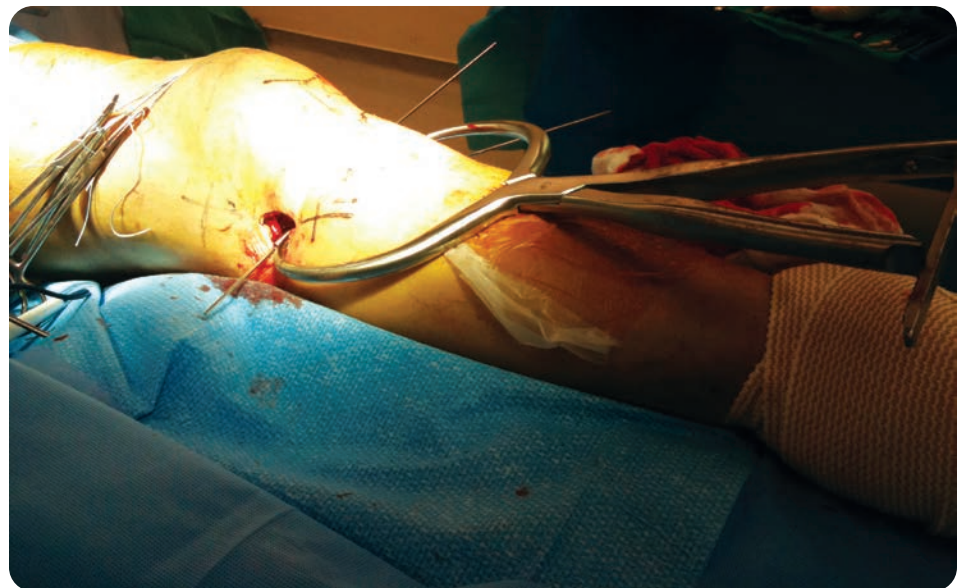

www.innomed.net

ISO 13485:2016

FREE TRIAL ON MOST INSTRUMENTS

INNOMED



103 Estus Drive, Savannah, GA 31404
info@innomed.net

912.236.0000 Phone
912.236.7766 Fax

Innomed-Europe Tel. +41 41 740 67 74
info@innomed-europe.com

1.800.548.2362

© 2023 Innomed, Inc.