

INNOMED

ORTHOPEDIC INSTRUMENTS



JULY
2026

Featuring many *New!* instruments inside



**McGrory
Foot Bolster**
Page 13



Arm Board
Page 17



**Cherf
Leg Holder**
Page 6



**Wixson Hip
Positioner**
Page 2

Patient Positioning

1.800.548.2362



INNOMED.NET

Wixson Hip Positioner

Designed by R.L. Wixson, MD

Provides stable positioning of a patient during hip surgery

The Wixson Hip Positioner is used for stable positioning of a patient during total hip and revision surgery. It is designed to be placed on top of the operating table.

The base plate is rubber-backed to reduce slipping on the table. The uprights can easily be slid in and out of the multiple slots in the plate for desired positioning and locked into position with the locking bolt. The complete upright assembly is radiolucent.

The upright pads and the base plate pad are made of semi-dense foam to help prevent pressure points and are sealed with a washable coating. The coating also helps to lessen the possibility of skin breakdown.

The hip positioner consists of: One 10" post with double pads, one 6" post with a single pad, one 20" base plate, one base plate pad, two 2" spacers, one 4" knob, and one 6" knob.

The spacers and longer knobs are supplied for use with larger patients: use one spacer with the 4" knob, or combine the two spacers to use with the 6" knob.

The pad assembly can be adjusted for additional height and width. The upright posts are modular. The complete unit is radiolucent, and all parts are autoclavable except for the washable foam pads.



Baseplate Dimensions: 20" x 11.25" (50,8 cm x 28,6 cm)

PRODUCT NO'S:

4050 [Set]
Consists of: One 10" post assembly with double pads, one 6" post assembly with a single pad, one 20" base plate, one base plate pad, two 2" spacers, one 4" knob, and one 6" knob.

Pads, Pad Parts, Pad Assemblies & Sets:

4050-LPD [Large Pad]
4150-CP [Curved Pad]
4150-P [Pad Plate] (Small Pad Attaches to this Plate)
4150-PA [Pad Assembly] (Pad Knob, Pad Plate & Pad)
4150-PB [Pad Knob] (Attaches Pad Plate to Post)
4150-PD3 [Set of 3 Small Pads]
4150-PD12 [Case of 12 Small Pads]
4150-V [Set of 30 Velcro Strips] (For Pad Attachment)

Posts, Post Parts, Post Assemblies & Sets:



4150-06 [6" Post (15,2 cm)] (Post Only)
4150-06B [6" Post with Pad Assembly]
4050-06C [6" Post with Pad Assembly and Slider]
4150-08 [8" Custom Post (20,3 cm)] (Post Only)
4150-09 [9" Custom Post (22,9 cm)] (Post Only)
4150-10 [10" Post (25,4 cm)] (Post Only)
4150-10B [10" Post with Pad Assembly]
4050-10C [10" Post with Pad Assembly and Slider]
4150-12 [12" Custom Post (30,5 cm)] (Post Only)
4150-14 [14" Custom Post (35,6 cm)] (Post Only)
4050-PA [Curved Pad Post Assembly] (Curved Pad, Pad, (2) Pad Knobs, (2) Pad Plates, 8" Post)
4050-PDS [Full Pad Set] (1 Large & 3 Small Pads)
4150-PS [Post Screw]
4050-S [Slider] (Connects Post to Baseplate)

Extension Parts, Extension Assemblies & Sets:


4150-C [2" Spacer (5,1 cm)]
4150-C4 [4" Spacer (10,2 cm)]
4150-EK [4" Knob (10,2 cm)]
For use with 2" Spacer
4150-EK4 [6" Extension Long Knob (15,2 cm)]
For use with two 2" Spacers or one 4" Spacer
4150-EK6 [8" Extension Long Knob (20,3 cm)]
For use with one 2" Spacer and one 4" Spacer
4150-EXT [2" Extension] (2" Spacer with 4" Knob)
4150-EXT4 [4" Extension] (4" Spacer with 6" Knob)
4150-EXT6 [6" Extension] (4" and 2" Spacers with 8" Knob)

Baseplates:


4050-BP [20" Standard Baseplate (50,8 cm)] (Baseplate Only)
4050-BP24 [24" Custom Baseplate (61 cm)] (Baseplate Only)

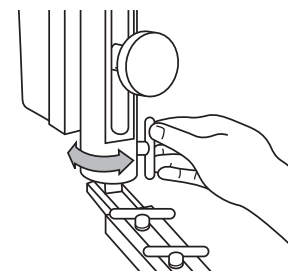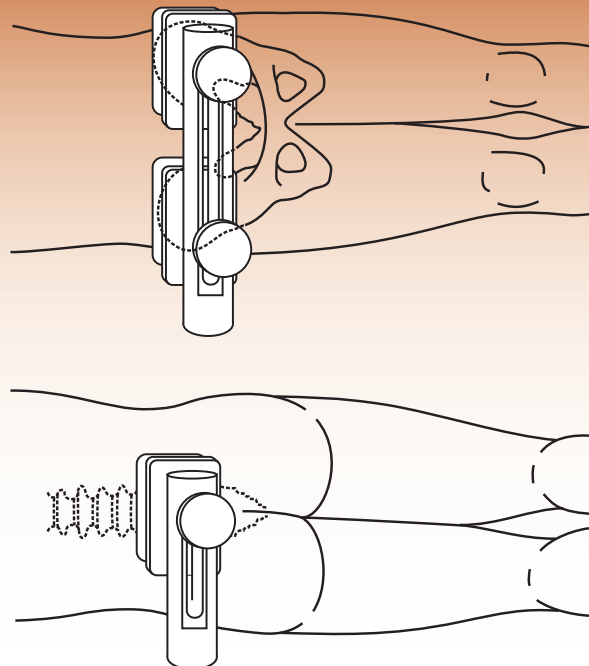

The upright assembly can be rotated and locked in place.

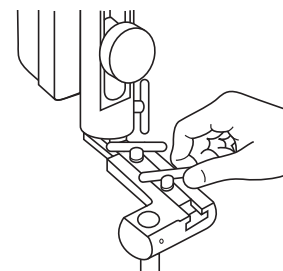

The upright assembly adjusts by a sliding track to accommodate various sized patients. It is locked in the sliding track by tightening one or two locking bolts.

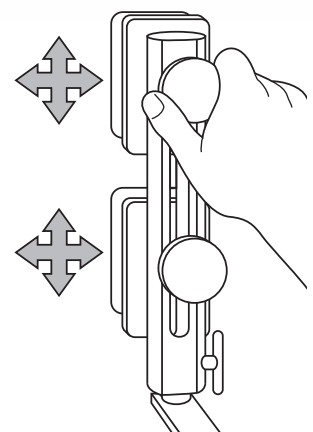

The pads can be adjusted for height and width.



Stulberg Hip Positioner

Designed by S. David Stulberg, MD

Provides stable positioning of a patient during hip surgery


The Stulberg Hip Positioner is used for stable positioning of a patient during total hip and revision surgery. It is designed to attach directly to the operating table utilizing the existing table adapters.

The upright pads are made of semi-dense foam to help prevent pressure points and are sealed with a washable coating. The coating also helps to lessen the possibility of skin breakdown.

The hip positioner consists of: One 10" post assembly with double pads and one 6" post assembly with a single pad, two 2" spacers, one 4" knob, one 6" knob, and two table attachments.

The spacers and longer knobs are supplied for use with larger patients: use one spacer with the 4" knob, or combine the two spacers to use with the 6" knob.

The pad assembly can be adjusted for additional height and width. The upright posts are modular. The complete unit is radiolucent, and all parts are autoclavable except for the washable foam pads and storage case.

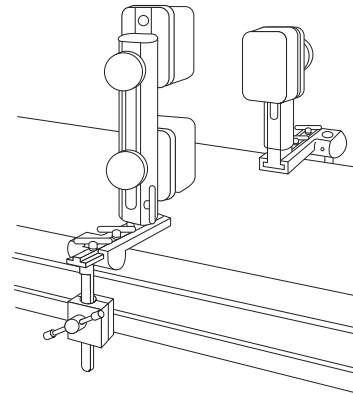



Storage Case Included

PRODUCT NO.'S:	
4150-00	[Set] Consists of: Case, one 10" post assembly with double pads, one 6" post assembly with a single pad, two 2" spacers, one 4" knob, one 6" knob, and two table attachments.
Pads, Pad Parts, Pad Assemblies & Sets:	
4150-CP	[Curved Pad]
4150-P	[Pad Plate] (Small Pad Attaches to this Plate)
4150-PA	[Pad Assembly] (Pad Knob, Pad Plate & Pad)
4150-PB	[Pad Knob] (Attaches Pad Plate to Post)
4150-PD3	[Set of 3 Small Pads]
4150-PD12	[Case of 12 Small Pads]
4150-V	[Set of 30 Velcro Strips] (For Pad Attachment)
Posts, Post Parts, Post Assemblies & Sets:	
4150-06	[6" Post (15,2 cm)] (Post Only)
4150-06B	[6" Post with Pad Assembly]
4150-06C	[6" Post with Pad Assembly and Slider]
4150-06D	[6" Post with Pad Assembly, Slider & Table Attachment]
4150-08	[8" Custom Post (20,3 cm)] (Post Only)
4150-09	[9" Custom Post (22,9 cm)] (Post Only)
4150-10	[10" Post (25,4 cm)] (Post Only)
4150-10B	[10" Post with Pad Assembly]
4150-10C	[10" Post with Pad Assembly and Slider]
4150-10D	[10" Post with Pad Assembly, Slider & Table Attachment]
4150-12	[12" Custom Post (30,5 cm)] (Post Only)
4150-14	[14" Custom Post (35,6 cm)] (Post Only)
4050-PA	[Curved Pad Post Assembly] (Curved Pad, Pad, (2) Pad Knobs, (2) Pad Plates, 8" Post)
4150-PS	[Post Screw]
4150-S	[Slider] (Connects Post to Table Attachment)
4150-TA	[Table Attachment]
Extension Parts, Extension Assemblies & Sets:	
4150-C	[2" Spacer (5,1 cm)]
4150-C4	[4" Spacer (10,2 cm)]
4150-EK	[4" Knob (10,2 cm)] For use with 2" Spacer
4150-EK4	[6" Extension Long Knob (15,2 cm)] For use with two 2" Spacers or one 4" Spacer
4150-EK6	[8" Extension Long Knob (20,3 cm)] For use with one 2" Spacer and one 4" Spacer
4150-EXT	[2" Extension (5,1 cm)] (2" Spacer with 4" Knob)
4150-EXT4	[4" Extension (10,2 cm)] (4" Spacer with 6" Knob)
4150-EXT6	[6" Extension (15,2 cm)] (4" and 2" Spacers with 8" Knob)
Other Parts:	
9002	[Storage Case]



Optional Hip Positioner Parts:

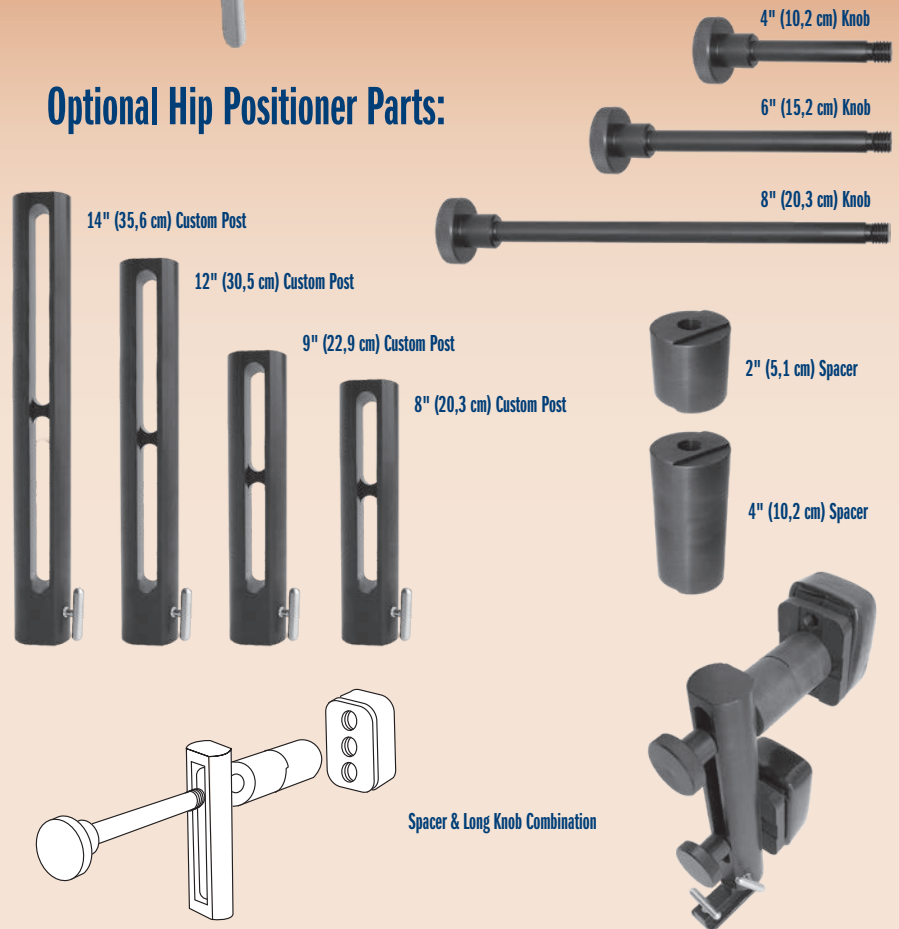


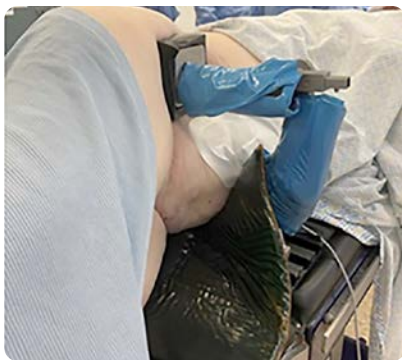


Belfast Sagittal Plane Positioner

Designed by David Beverland, FRCS

A sturdy and stable patient support system for posterior approach total hip arthroplasty in the lateral decubitus position

- ▶ Does not attach to the table, making it compatible with all OR tables
- ▶ Very secure and easy to tighten
- ▶ Accommodates the very obese patient

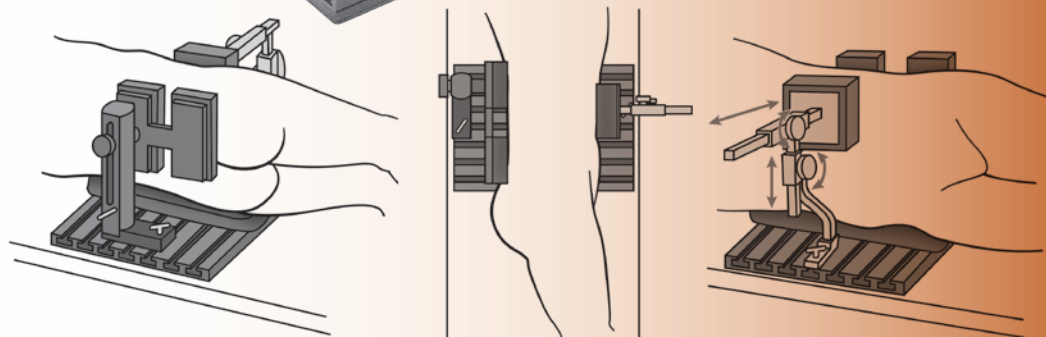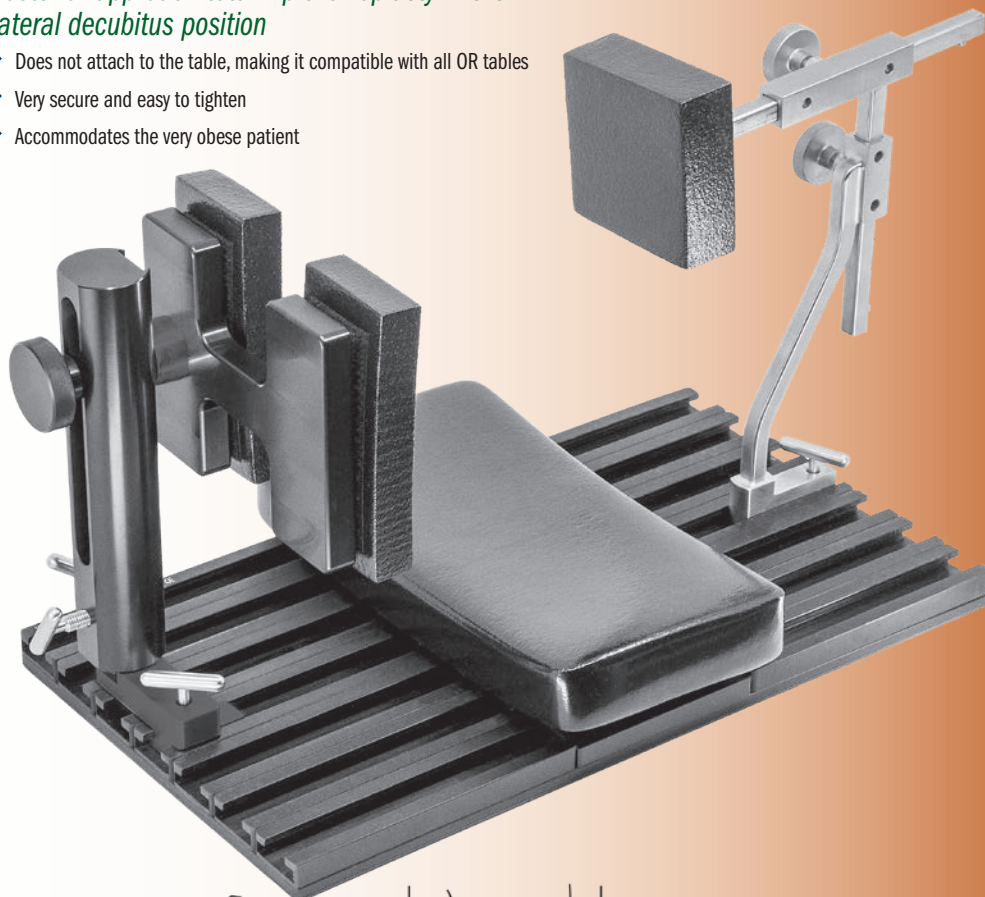

- ▶ Anterior ASIS support unit allows both height and depth adjustment, eliminating the need to adjust the baseplate locking mechanism once secured. The depth adjustment allows extension of up to 150 mm beyond the support base, providing space to avoid the abdominal apron.
- ▶ Posterior support unit offers both Sacral and Iliac Crest support, connected together as one unit, and provides both height and rotation adjustments, as well as a hinged pivot where the iliac crest side attaches the unit to the post, that can be loosened and adjusted to allow rotation in 5° increments.
- ▶ Allows access to the PSIS before and during surgery. This can help to ensure that the sagittal plane of the pelvis is horizontal at the time of cup implantation.

PRODUCT NO'S:	
4170-00	[Complete Set]
Set Includes / Available Individually:	
4170-03	[Anterior Upright Support]
4170-04	[Anterior Plane Support]
4170-05	[Anterior Clamp Support]
4170-AKS	[Anterior Knob Screw] <i>Two (2) included in Set, One (1) with this product number</i>
4170-AP	[Anterior Plane Pad]
4170-06	[Posterior Sagittal Plane Support]
4170-07	[Posterior 9.5" Post]
4170-PKS	[Posterior Knob Screw]
4170-08	[Posterior Base]
4170-09	[Posterior Angle Adjuster]
4170-T	[Posterior T-Handle Screw]
4170-PP	[Posterior Support Pad] <i>Two (2) included in Set, One (1) with this product number</i>
4150-PS	[Post Screw] <i>Three (3) included in Set, One (1) with this product number</i>
4050-BP	[20" Baseplate Only]
4050-LPD	[Hip Positioner Large Pad]




Extension Set for Anterior THR Tables

Designed by David Ott, MD

Designed to add lift to the femoral hook during an anterior THR case and be able to remove without breaking the sterile field



PRODUCT NO'S:	
8004-00	[Set of One Each]
Also available individually:	
8004-S	[Short Extension] Extension Length: 2" (5,1 cm) Overall Length: 2.6" (6,6 cm)
8004-L	[Long Extension] Extension Length: 3" (7,7 cm) Overall Length: 3.625" (9,2 cm)

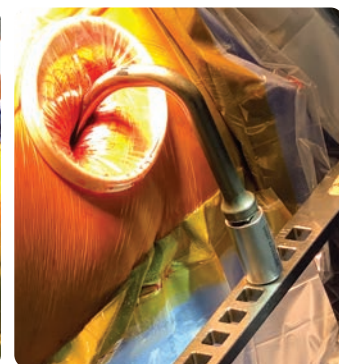



Thornberry Hip Positioner

Designed by Robert L. Thornberry, MD

Designed to be adjustable yet sturdy, and is especially helpful when stabilizing a large patient during total hip and revision surgery

The Thornberry Hip Positioner is designed to attach directly to the operating table utilizing existing table clamps, or the Innomed #2595 Table Clamps, which are not included.

The upper arm assembly can be adjusted for height. Both arms include a push-button to allow the pad platform to swivel and lock into any of three fixed positions. The tall 18" post also includes a push-button to allow the post/arms unit to swivel and lock into any of three fixed positions.

The complete unit is autoclavable except for the foam pads. The pads are made of semi-dense foam to help prevent pressure points and are sealed with a washable coating. The coating also helps to lessen the possibility of skin breakdown.

It may be necessary to place the Double Table Attachment(s) 180°—sticking out from the table—to accommodate the large patient, as shown at left.

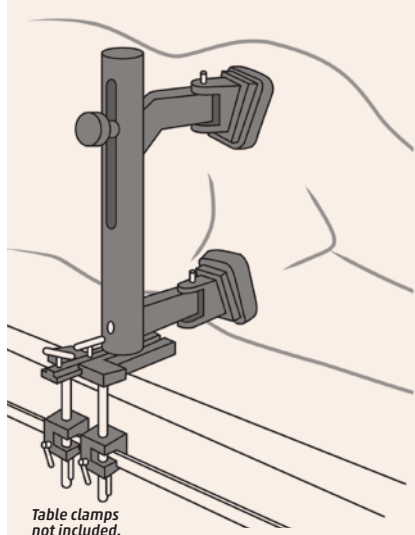


Table clamps not included.


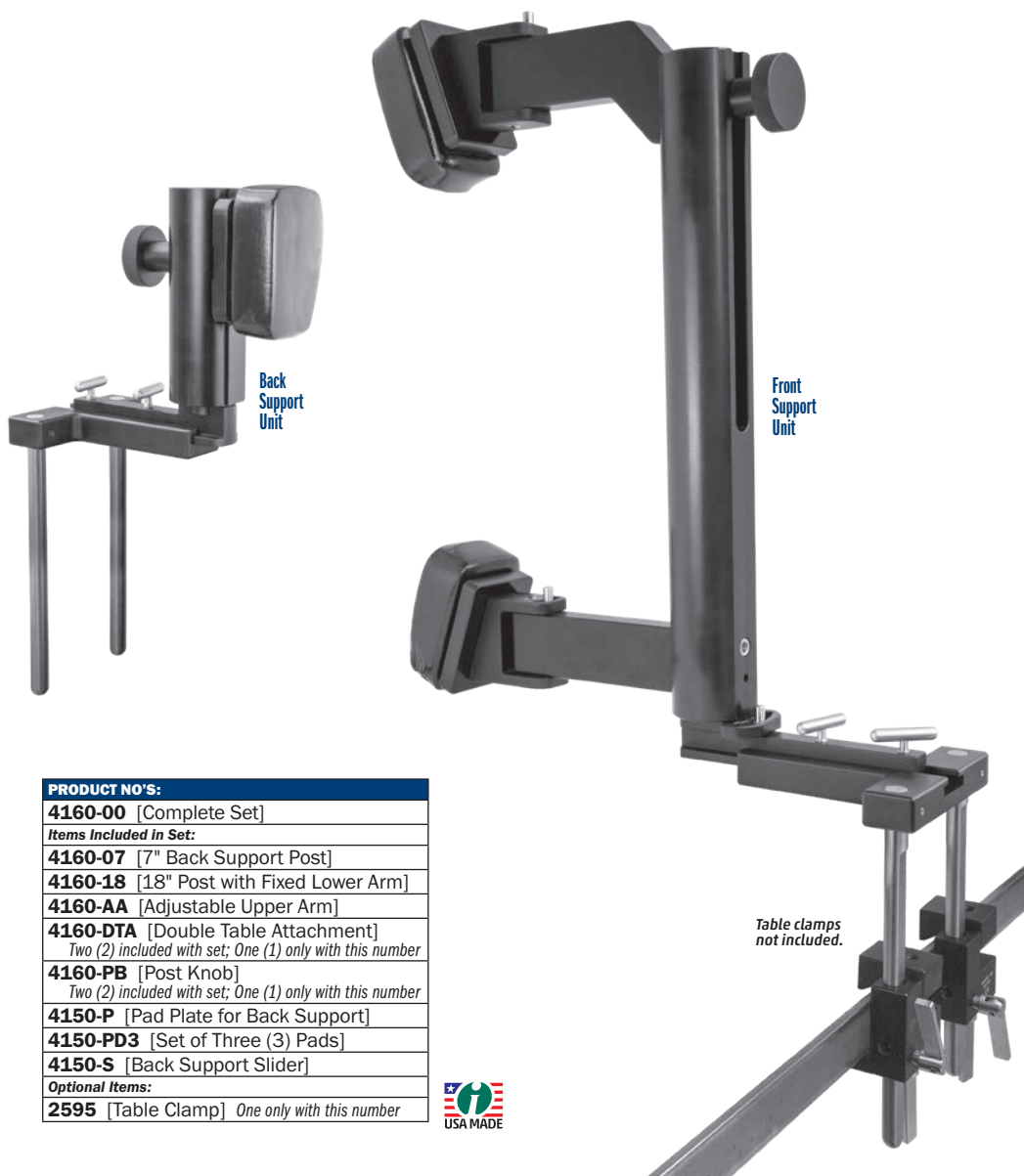


Table clamps not included.

PRODUCT NO.'S:	
4160-00	[Complete Set]
Items Included in Set:	
4160-07	[7" Back Support Post]
4160-18	[18" Post with Fixed Lower Arm]
4160-AA	[Adjustable Upper Arm]
4160-DTA	[Double Table Attachment]
Two (2) included with set; One (1) only with this number	
4160-PB	[Post Knob]
Two (2) included with set; One (1) only with this number	
4150-P	[Pad Plate for Back Support]
4150-PD3	[Set of Three (3) Pads]
4150-S	[Back Support Slider]
Optional Items:	
2595	[Table Clamp] One only with this number

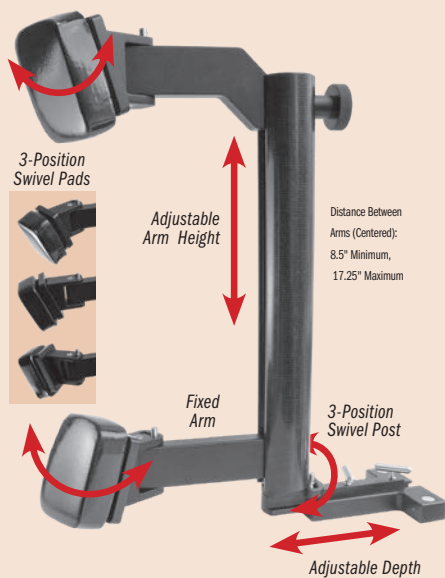

Positioner consists of:

Front Support Unit

One 18" post assembly with a lower arm and swiveling pad platform, one upper arm assembly with a swiveling pad platform, one post knob, two pads, and one double table attachment.

Back Support Unit

One 7" post assembly, one post knob, one pad plate, one pad, one slider, and one double table attachment.



- ▶ Allows one assistant to secure the leg for femoral preparation
- ▶ Attaches directly to a standard operating table
- ▶ Allows easy assessment of hip stability and leg length discrepancy
- ▶ Scale on the support rod helps to allow for precise and reproducible placement of the leg positoner according to surgeon preference

See page 19 for Reference Measuring Statement for devices with scales, markings, or visual indicators.

US Patent No. 11,744,757




Direct Anterior Total Hip Arthroplasty Leg Positoner Assembly

Designed by Benjamin M. Frye, MD

Designed to help position the operative leg for direct anterior approach total hip arthroplasty using a standard operating table

PRODUCT NO'S:	
4165-00	[Assembly with Standard Table Clamp]
4165-00-EU	[Assembly with EU Table Clamp]
Individual/Replacement Parts:	
4165-01	[Leg Support] Overall Length: 8.75" (22,2 cm) Width: 2.5" (6,4 cm)
4165-02	[Support Rod] Dimensions: 9.625" x 7.125" x 1" (24,5 x 18,1 x 2,5 cm)
4165-P	[Leg Support Pad] Dimensions: 8.5" x 3" x 1.75" (21,6 x 7,6 x 4,4 cm)
4150-PS	[Support Screw] T-Handle Width: 1.875" (4,7 cm) T-Handle Depth: 1.5" (3,8 cm)
9120	[Table Clamp]
9120-EU	[Table Clamp for European Bed Rails]

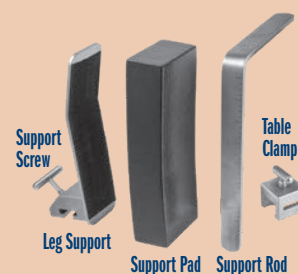


Das Anterior Hip Bolster Assembly

Designed by Amal Das, MD

Designed to help provide counter resistance on the contralateral hip during reaming and implant insertion in direct anterior hip arthroplasty

PRODUCT NO'S:	
4166-00	[Das Anterior Hip Bolster Assembly]
Set Includes / Available Individually:	
4166-01	[Das Anterior Hip Bolster Support]
4166-02	[Das Anterior Hip Bolster Rod]
Set Includes / Replacement Parts:	
2595	[Table Clamp]
4150-PD2	[Positioning Pads - Set of 2]
4150-PS	[Post Screw]

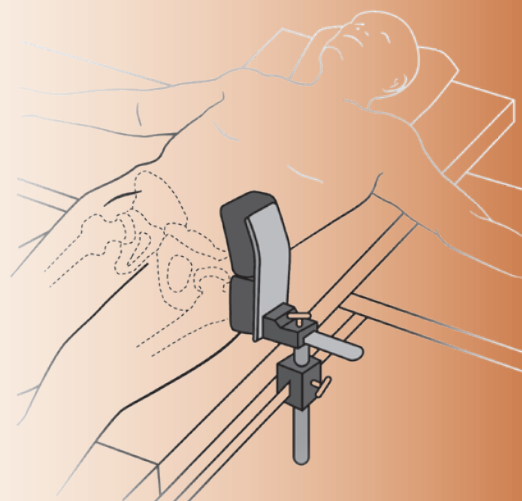

Cherf Leg Holder

Designed by John Cherf, MD

Supports the lower extremity for prepping before hip or knee surgery

Useful for all lower extremity procedures and is particularly helpful for supporting the leg with the patient positioned in the lateral position. By holding the foot/ankle in an externally rotated position, the knee can be locked into extension which helps eliminate the need for manual support.

The leg holder may also be used to support the limb for surgical patients in the supine position such as for knee and foot/ankle procedures.

PRODUCT NO'S:	
2270	
Replacement Parts:	
4150-PD3	[Set of 3 Small Pads]



Capello Patient Positioner

Designed by William Capello, MD

Provides stable positioning of a patient during hip procedures

Two-piece board design available for easy handling
Interlocking board pieces



Also available in a one-piece design



PRODUCT NO'S:

4090 [Set with 2-Piece Board and Standard Table Clamp]

4090-EU [Set with 2-Piece Board and EU Table Clamp]

4095 [Set with 1-Piece Board and Standard Table Clamp]

4095-EU [Set with 1-Piece Board and EU Table Clamp]

Replacement Parts:

4090-PB [2-Piece Positioning Board]

4095-PB [1-Piece Positioning Board]

4090-06 [6" (15,2 cm) Radiolucent Peg]

Four included in set; one with this product number

4090-08 [9" (22,9 cm) Radiolucent Peg]

Four included in set; one with this product number

4090-SC [Stabilizing Clamp]

Two included in set; one with this product number

4090-01 [Large Gel Pad]

9120 [Table Clamp]

Two included in set; one with this product number

9120-EU [Table Clamp for European Bed Rails]

Two included in set; one with this product number

Optional Parts:

4090-02 [Peg Gel Pad]

4090-EXT2 [Peg Extension - 2" (5,1 cm)]

4090-EXT [Peg Extension - 4" (10,2 cm)]

4090-EXT6 [Peg Extension - 6" (15,2 cm)]

4090-EXT8 [Peg Extension - 8" (20,3 cm)]

6" (15,2 cm) & 9" (22,9 cm) Radiolucent Pegs

Gel Pad

Stabilizing Clamp

Table Clamp

Board

Board Dimensions:
47" x 18.75" (119,4 cm x 47,6 cm)

Each set includes:

Board, Gel Pad,
(4) 6" Radiolucent Pegs,
(4) 9" Radiolucent Pegs,
(2) Stabilizing Clamps,
(2) Table Clamps.


Optional two-piece board construction allows for easier use and storage. The board is also available in a one-piece design. All gel pads, pegs and peg height extensions can be used with existing peg boards.



Optional Peg Pad



Optional
2", 4", 6", & 8"
Peg Extensions

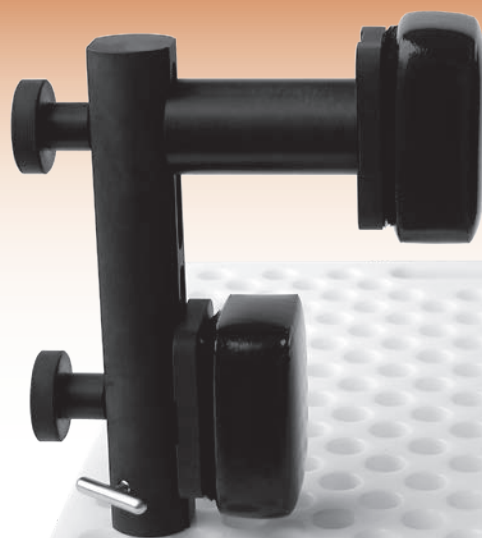

Large Patient Peg Board Positioner Post Assembly

Designed by Paul Ramsey, MD

Especially helpful with large patients where reaching the a.s.i.s. is needed for stabilization

PRODUCT NO'S:

4150-10P [Complete Set]

Assembly Set Includes:

4090-03 [Post Assembly Adapter]

4150-10B [10" (25,4 cm) Post with 2 Pads]

4150-EXT [2" (5,1 cm) Spacer with 4" (10,2 cm) Knob]

4150-EXT4 [4" (10,2 cm) Spacer with 6" (15,2 cm) Knob]

NOTE: The peg board positioner is available separately and is not included with this assembly set.



Post Assembly Adapter





Multi-Adjustment Hip Positioner

Provides stable positioning of a patient during hip surgery

Multi-adjustment arms allow the positioner to be adjusted to fit all sizes of patients. Extra attachment allows for more versatility of placement. Especially helpful with large patients where reaching the a.s.i.s. is needed for stabilization.

PRODUCT NO'S:

4030

Replacement Part:

4150-PD2 [Set of 2 Small Pads]



Use in Knee Surgery

Use in Hip Surgery

Berger Block Positioner

Designed by Richard Berger, MD

Designed for lower extremity positioning with dual height options



PRODUCT NO'S:

2750-00 [Set]

Dimensions with Pads: 4.75" x 6.75" x 8"
(12,1 x 17,1 x 20,3 cm)

Set Includes / Available Individually:

2750-01 [Positioner Only]


Dimensions: 4.125" x 6.125" x 8"
(10,5 x 15,6 x 20,3 cm)

2750-P [Positioner Pad Only]

[Brown Strap Only] (2) Included in Set

Optional Items:

2750-S [Brown Straps] **Pkg of 10**

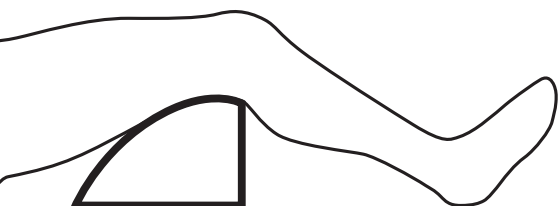

Lower Extremity Leg Positioner

Designed by Ronald Romanelli, MD

Used to support knee and leg during surgery, and can be used for casting

- ▶ Utilized for rodding of femurs or tibias
- ▶ Also useful for knee surgery and closures
- ▶ Very supportive, distributes stresses on leg, used instead of bolsters
- ▶ Supplied with one autoclavable silicone pad
- ▶ Aluminum positioner is radiolucent and gas or steam sterilizable

PRODUCT NO'S:

2745

Dimensions: 5.5" H x 9.5" L x 9.25" W
(12,7 cm x 24,1 cm x 23,5 cm)

Replacement Parts:

2760-P [Silicone Pad]




Stulberg Sliding Bolster

Designed by S. David Stulberg, MD

Helps eliminate the need for a sand bag during total knee surgery

The base plate is attached to the table and the sterile sliding bolster is placed on top of the sterile drape. The bolster can be adjusted for different angles of knee flexion during surgery.

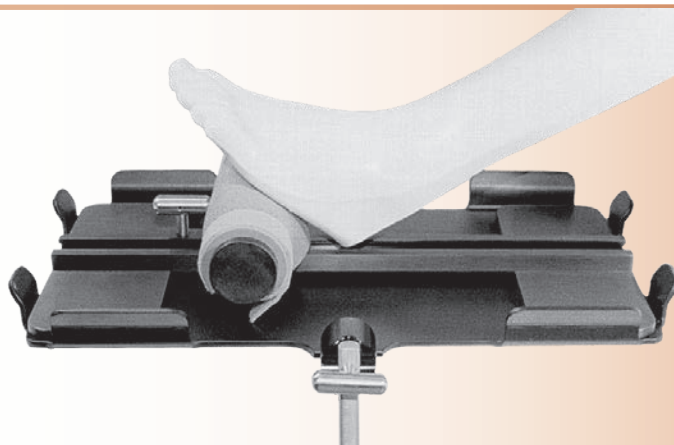


PRODUCT NO:

2730

Base Dimensions:
20" x 10.5" (50,8 cm x 26,7 cm)

Patient Self Stress Assembly Set

Designed by Kyle Cook, RTR and David Mauerhan, MD

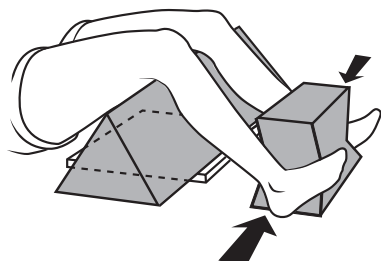

Designed to help position a patient for X-ray evaluation to help determine candidacy for Unicondylar Knee Arthroplasty

Triangle Positioner
#2741-01



Contoured Cube
#2741-02




PRODUCT NO'S:	
2741-00	[Set]
Individual Instruments:	
2741-01	[Triangle Positioner] Dimensions: 24" x 9" x 9" (61 cm x 23 cm x 23 cm)
2741-02	[Contoured Cube] Dimensions: 11" x 9" x 6" (28 cm x 23 cm x 15,2 cm)



Chandran Thigh Lift Positioner

Designed by Rama E. Chandran, MD

Designed to help lift and position the thigh from above during knee surgery

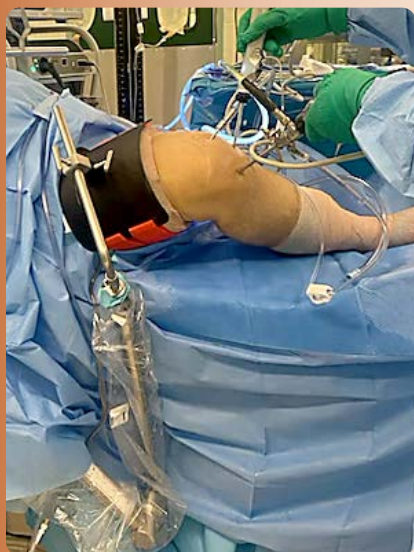

The optional thigh lift adapter is designed for use with a hydraulic lift device instead of the manual lift rod with table clamp.



PRODUCT NO'S:	
4167-00	[Positioner Set]
Set Includes:	
4167-01	[Cradle] Dimensions: 7.5" x 7.75" (19,1 x 19,7 cm)
4167-02	[Rod] Width: 13.6" (34,5 cm) Height: 20.6" (52,3 cm)
4150-PS	[Post Screw] Two included in set, one with this product number
2760-P	[Silicone Pad] Dimensions: 9" x 8" x 1/2" thick (22,9 x 20,3 x 1,3 cm)
[Orange Strap] Dimensions: 2" x 31.5" (5,1 x 80 cm) Two included in set, see below for replacement 10-pack	
[Short Black Straps] Dimensions (Each): 2" x 24" (5,1 x 61 cm) Two included in set, see below for replacement 10-pack	

PRODUCT NO'S:	
Replacement Parts:	
2760-P	[Silicone Pad] Dimensions: 9" x 8" x 1/2" thick (22,9 x 20,3 x 1,3 cm)
8129-P	[Short Black Straps] 10-Pack Dimensions (Each): 2" x 24" (5,1 x 61 cm)
8130-P	[Orange Straps] 10-Pack Dimensions (Each): 2" x 31.5" (5,1 x 80 cm)

PRODUCT NO'S:	
Optional Parts:	
4167-03	[Adapter] Width: 12.4" (31,5 cm) Height: 4.875" (12,4 cm)
9125	[Rotating Table Clamp]



Durham Leg Positioner

Designed by Al Durham, MD

Placed against the thigh, helping to hold the leg upright in knee surgery

Supplied with a sterilizable table clamp. The pad is made of semi-dense foam to help prevent pressure points and is sealed with a washable coating.

PRODUCT NO'S:

4105 [Set with Standard Table Clamp]


4105-EU [Set with EU Table Clamp]

Replacement Parts:


4105-P [Pad]

9120 [Table Clamp]

9120-EU [Table Clamp for European Bed Rails]

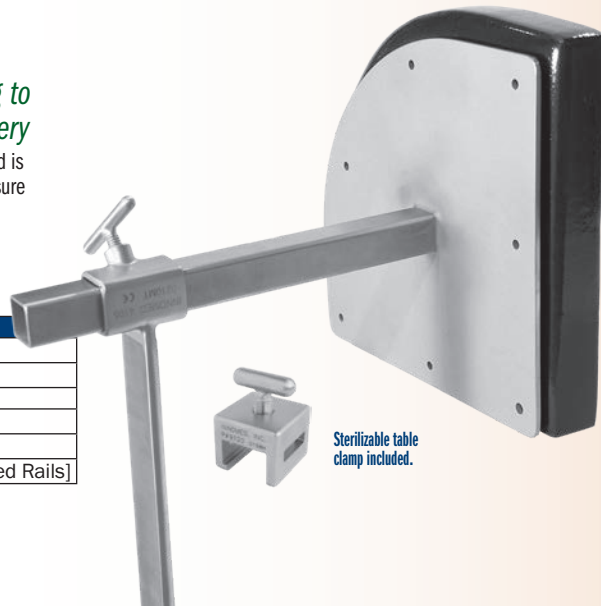

Sterilizable table clamp included.

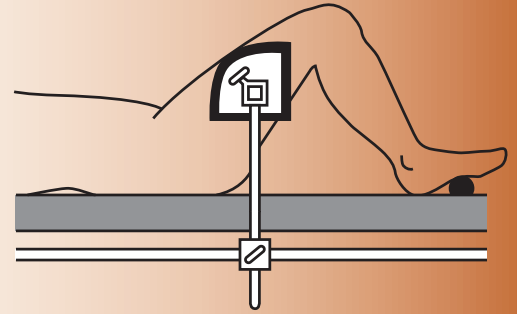

Modified 90° Leg Stabilizer

Designed by Gregory Fanelli, MD

Useful in total knee surgery to hold the leg in position

Helps to open up the knee joint when pressure is applied to the lower leg. Pad and sterilizable table clamp included.

PRODUCT NO'S:

2725 [Set with Standard Table Clamp]

Post Height: 11.375" (28.9 cm)

Pad Length: 9" (22.9 cm)

Pad Diameter: 3" (7.6 cm)


2725-EU [Set with EU Table Clamp]

Replacement Parts:


8840-P [Pad]

9120 [Table Clamp]

9120-EU [Table Clamp for European Bed Rails]

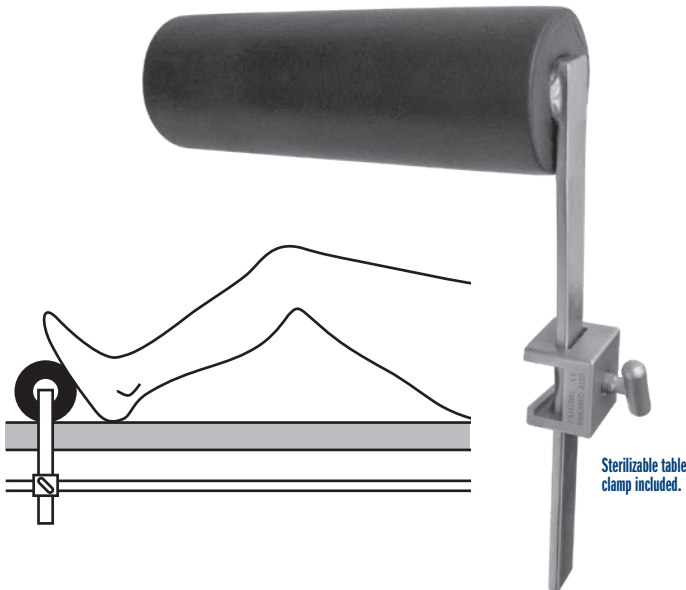


Sterilizable table clamp included.

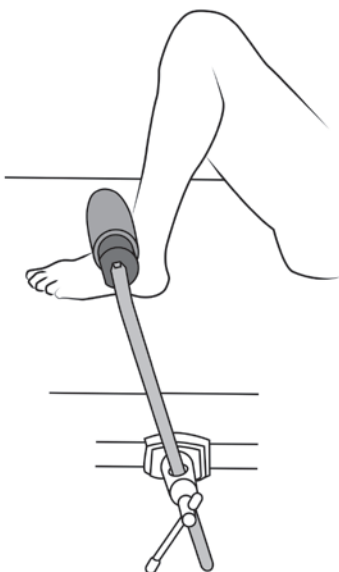

Hyperflex Foot Positioner Assembly

Designed by Morteza Meftah, MD and Ira Kirschenbaum, MD

Designed to help secure the foot for positioning of the knee in the hyperflex position

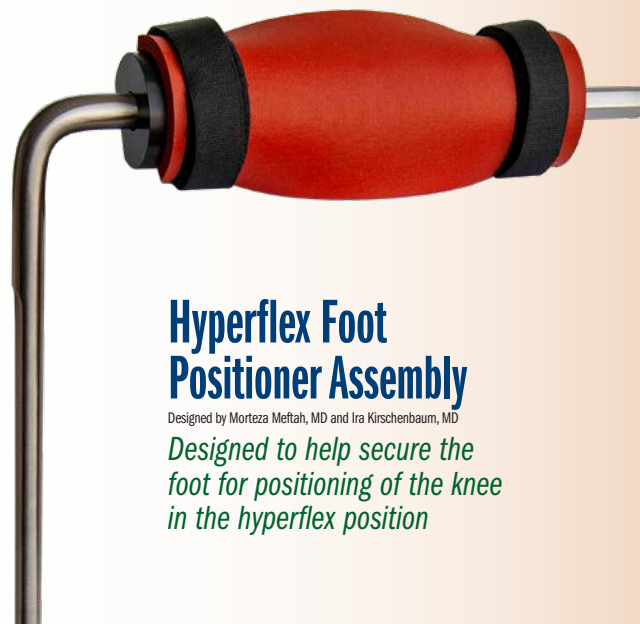

PRODUCT NO'S:

2589-00 [Complete Assembly]
Overall Length: 19" (48,3 cm)

Individual/Replacement Parts:

2589-01 [Foot Positioner]
Overall Length: 18" (45,7 cm)

2590-B [Positioner Bar]
Overall Length: 19" (48,3 cm)

2730-P [Pad & Two Straps]

4150-PS [Post Screw]

Optional Part:

2590-S [Black Straps] Pkg of 10



FREE TRIAL

Most instruments a no-charge two-week trial includes FREE UPS

*When shipped to a hospital, an additional charge applies. Free trial offer excludes instruments, which are available for purchase at a replacement charge.

Stanton Arthroscopic Leg Holder

Designed by John Stanton, MD

- ▶ Sliding leg holder can be adjusted for small calves or to accommodate large thighs
- ▶ Locking pin prevents sides from spreading apart
- ▶ Strap can be placed high or low through the slots in the side plates to accommodate large/small limbs
- ▶ Strap is strongly secured with a toothed clamp
- ▶ Support rod, when clamped into a standard table clamp, helps to prevent rotation

PRODUCT NO'S:


4045

Dimensions: 16.5" L x 8.5" H x 3.5" W
(42 cm L x 21,6 cm H x 9 cm W)
Fits Legs: From 4" to 11" (10 cm to 28 cm)

Replacement Parts

4045-S [Strap]


Overall Length: 28" (71,2 cm)

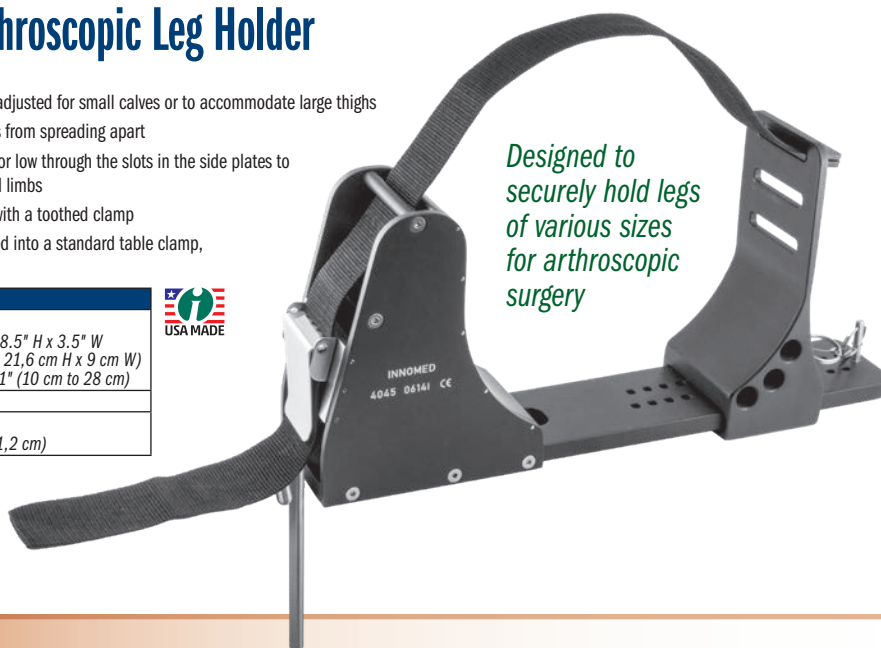

Designed to securely hold legs of various sizes for arthroscopic surgery



PROTOTYPE SHOWN

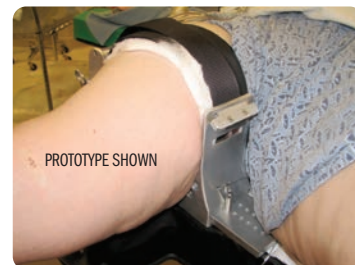


PROTOTYPE SHOWN



are available for
week evaluation –
Ground Shipping*

hospital or medical center;
s for expedited shipping.
ies implant extraction
available as rentals. There is a
with the hip positioners.

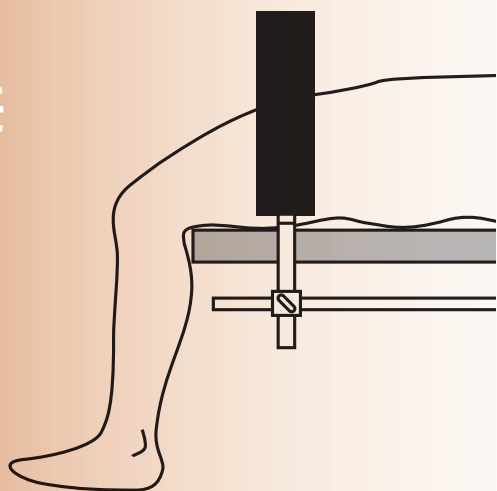

Leg Stabilizer

Designed by Gregory Fanelli, MD

Useful in arthroscopic knee surgery to hold the leg in position

Helps to open up the knee joint when pressure is applied to the lower leg. Sterilizable table clamp included.

PRODUCT NO'S:

8840 [Set with Standard Table Clamp]

Overall Length: 18.5" (47 cm)
Handle Length: 9.25" (23,5 cm)
Pad Diameter: 3" (7,6 cm)

8840-EU [Set with EU Table Clamp]

Replacement Parts:

8840-P [Pad]

9120 [Table Clamp]

9120-EU [Table Clamp for European Bed Rails]



Sterilizable table
clamp included.



George Arthroscopic Knee Positioner

Designed by Michael S. George, MD

Provides lateral as well as superior support which allows valgus stress to open the medial compartment

Shape does not squeeze the thigh, making the need for a thigh tourniquet optional. If desired, the unit can easily be rotated out of the way without disrupting the sterile field. Using with a standard operating table clamp, the unit can easily be raised or lowered to accommodate all thigh sizes.


PRODUCT NO'S:

2735

Overall Height: 22" (55,9 cm)
Post Height: 12" (30,5 cm)
Pad Width: 3" (7,6 cm)

Replacement Part:

2735-P [Pad]

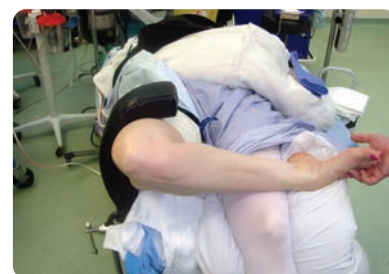



Adjustable Knee & Tibial Positioner

Designed by Ashutosh Chaudhari, MD

Adjustable design allows for use in procedures around the knee such as tibial nailing, tibial condyle plating, patella fracture fixation, supracondylar fracture plating, supracondylar fracture nailing, and total knee replacement

PRODUCT NO'S:

2770-00 [Set]
Includes Positioner, Pad,
and Two Short Straps

Individual / Replacement Parts

2770-01 [Positioner]
Overall Length (Folded): 28" (71,1 cm)
Overall Length (Flat): 54.75" (139 cm)
Maximum Triangle Height: 14" (35,6 cm)
Width: 5.5" (14 cm)
Thickness (Folded): 1.8" (4,6 cm)
Thickness (Flat): .75" (1,9 cm)

2770-P [Silicone Pad]
Dimensions: 12" x 5.5" (30,5 cm X 14 cm)

2590-S [Short Straps] Pkg of 10



Radiolucent. Steam sterilizable.



Sanders Extremity Positioning Tubes

Designed by Richard A. Sanders, MD

Designed to support the knee and ankle during lower extremity surgery

The 6" tube lifts the knee off the operating table and allows for approximately 30° of knee flexion. Very useful for closure of total knee incisions, supporting fractures of the distal femur, and tibia plateau fractures. The 4" tube elevates the foot and ankle for ankle fracture surgery. The tubes are made of aluminum, allowing them to be autoclaved. They help eliminate the need for rolled sheet bolsters.

PRODUCT NO'S:

2740-01 [Small]
Diameter: 4" (10,2 cm)
Width: 8" (20,3 cm)


2740-02 [Large]
Diameter: 6" (15,2 cm)
Width: 8" (20,3 cm)



Kirschenbaum Foot Positioner

Designed by Ira Kirschenbaum, MD

Helps eliminate the use of sand bags under the drape during total knee surgery

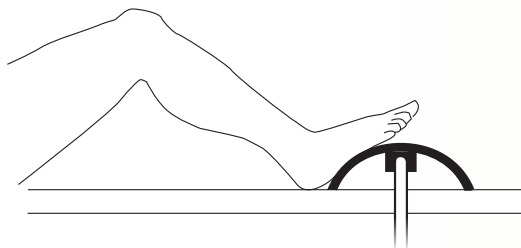

The foot rest is dome shaped for optimal foot contact and positioning the leg in flexion, and can be rotated. The unit can be used under the drape by attaching it to a standard table attachment or it can be sterilized for use on top of the drape. It can be attached to the table with the optional sterilizable table clamp. Supplied with a removable, sterilizable silicone foot pad.



Long #2590

Short #2591



PRODUCT NO'S:

2590 [Foot Positioner - Long] 15.5" x 6" (39,4 cm x 15,2 cm)

2591 [Foot Positioner - Short] 9.5" x 6" (24,1 cm x 15,2 cm)

Optional & Replacement Parts:

2590-P [Large Replacement Pad] 16" x 9" (40,7 cm x 22,9 cm)

2591-P [Small Replacement Pad] 9.5" x 9.25" (24,1 cm x 23,5 cm)

2595 [Optional Table Clamp]

Fromm Femur & Tibia Triangles

Designed by S.E. Fromm, MD *
Extra Small Triangle designed by S.E. Fromm, MD & Kenneth Merriman, MD

Used for femur and tibia positioning during nailing, repairs and fractures

Designed to position and hold the femur and tibia during intramedullary nailing of the tibia, ligament repairs and extremity fractures. Allows knee to be flexed greater than 90° to allow reaming and nail insertion without displacing fracture. The triangles are available in four heights: 8.5", 11", 14", and 16". The three smaller triangles are designed to fit inside the larger triangle for storage. They are supplied with an autoclavable silicone cushioning pad and velcro* straps. The triangles are also radiolucent and gas or steam sterilizable.

PRODUCT NO'S:	
2760-00	[Set of 3] Angles: Top 30°, Two Bottom 75°
Set Includes / Available Individually:	
2760-01	[11"] Base: 6" (15,2 cm), Height: 11" (27,9 cm)
2760-02	[14"] Base: 7" (17,8 cm), Height: 14" (35,6 cm)
2760-03	[16"] Base: 9" (22,9 cm), Height: 16" (40,7 cm)
Sold Separately - Not In Set:	
2760-XS	[8.5"] Base 5" (12,7 cm), Height: 8.5" (21,6 cm)
Replacement Parts:	
2760-P	[Silicone Pad]
2760-S	[Straps] Package of 18
8120-SP	[Straps for XS] Package of 10



*Velcro is a registered trademark of the Velcro Companies.




16" #2760-03

14" #2760-02


11" #2760-01

8.5" #2760-XS Sold Separately - Not In Set

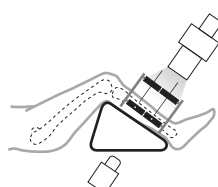


Tibia Reduced For:

- Open Reduction and Internal Fixation (ORIF)
- Application of uni- or multi-plane external fixator
- Knee ligament repairs and/or reconstruction

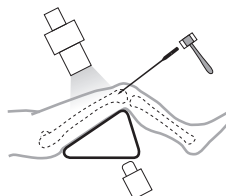

Retrograde Femoral Nailing

Triangle holds femur reduced (prevents sagging)

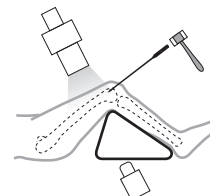

Retrograde Femoral Nailing



Tibial Nailing



McGrory Foot Bolster

Designed by James McGrory, MD

A foot bolster designed for knee surgery that allows for quick and easy adjustment of the flexion angle as surgery progresses



PRODUCT NO'S:	
2616-00	[McGrory Foot Bolster Assembly]
Assembly Includes / Available Individually:	
2616-01	[McGrory Foot Bolster] Overall Length: 10.65" (27 cm) Bolster Block Length: 9.5" (24,1 cm) Bolster Block Height: 3" (7,6 cm) Bolster Block Depth: 2.25" (5,7 cm)
2616-02	[McGrory Baseplate] Overall Length: 16.125" (41 cm) Overall Width: 10.875" (27,6 cm)
2620-01	[Tightening Screw] Overall Length: 2.25" (5,7 cm) Diameter: (16 mm)
2620-ARM	[Table Attachment] Overall Length: 12" (30,5 cm) Diameter: (16 mm)
2595	[Table Clamp]



Table Clamp Included in Assembly

- ▶ Easy adjustment of foot bolster position during surgery to place the knee at the desired flexion angle. Helpful as increasing flexion is gained as the surgery progresses.
- ▶ "Gas pedal" design of the foot bolster with anti-slide ridges to keep foot stable
- ▶ Bottom of the bolster pedal can be used to capture the great toe during maximum flexion and external rotation during difficult exposures and revisions
- ▶ Useful for holding the knee stable for robotic use

McGrory Foot Bolster


McGrory Baseplate

Connelly Leg Positioning Set

Designed by Camille Connelly, MD

A clamp free leg positioning board for partial or total knee replacement designed to go UNDER the drapes, with a sterilizable footpiece for easy positioning

New!

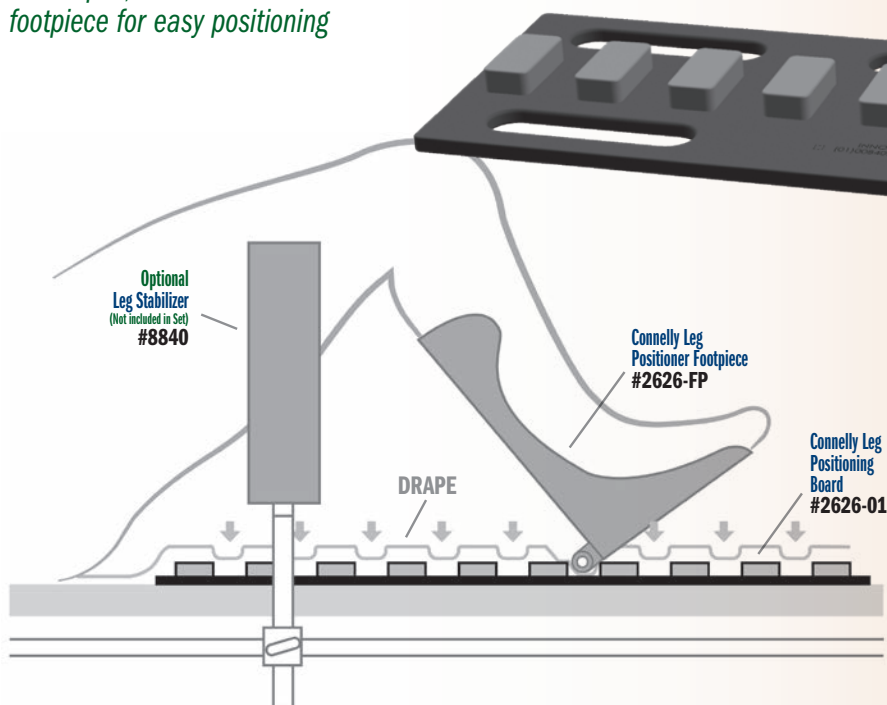

- ▶ The board can be stabilized on the OR table with Velcro or tape
- ▶ Two (2) Sterile Pads/Wraps are included with each new purchase (See right for Pad/Wrap information)
- ▶ Set includes the footpiece

PRODUCT NO'S:

2626-00 [Connelly Leg Positioning Set]

Set Includes / Available Individually:

2626-01 [Connelly Leg Positioning Board]

Overall Length: 30" (76,2 cm)

Overall Width: 10" (25,4 cm)

Height: 1.25" (3,2 cm)

2626-FP [Connelly Leg Positioner Footpiece]

Length: 11.9" (30 cm)

Height: 7.9" (20 cm)

Width: 6" (15,2 cm)

2629-L [Sterile Pad/Wrap]

Two included in set; product number is for one only

Replacement Parts:

2629-00 [Case of 10 Sterile Pads/Wraps]




Robb Leg Positioner

Designed by William Robb, MD

Provides stable positioning of the knee during surgery

- ▶ Slotted base allows the leg to be easily flexed or extended during knee surgery
- ▶ The complete unit is steam and gas sterilizable
- ▶ Supplied with a sterilizable table clamp which can be clamped over the sterile drape to the O.R. table side bar
- ▶ Three (3) Sterile Pads/Wraps are included with each new purchase (See right for Pad/Wrap information)

Slots are also designed to allow the foot piece to be rotated

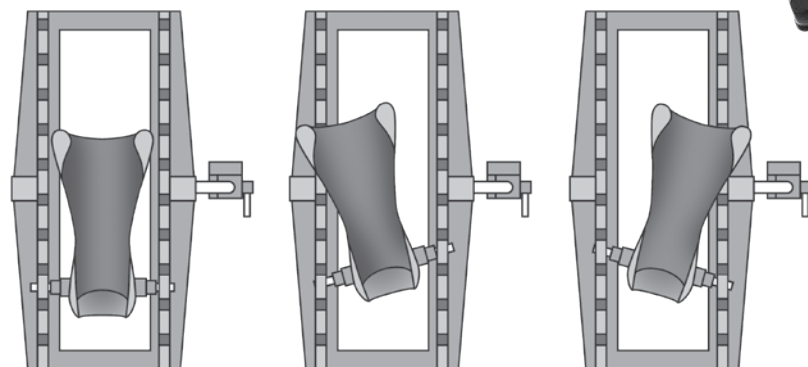

PRODUCT NO'S:

Base Dimensions: 21" x 11" (53,4cm x 27,9cm)

2630-11 [Robb Leg Positioner]

Optional & Replacement Parts:

2630-FP [Aluminum Footpiece Only]

2629-00 [Case of 10 Sterile Pads/Wraps]

2595 [Table Clamp]



Stulberg Leg Positioner


Designed by S. David Stulberg, MD

Provides stable positioning of the knee during surgery

PRODUCT NO'S:	
Base Dimensions: 20" x 10.5" (50,8 cm x 26,7 cm)	
2620-10	[Stulberg Leg Positioner]
Optional & Replacement Parts:	
2620-FP	[Aluminum Footpiece Only]
2629-00	[Case of 10 Sterile Pads/Wraps]




Usage guide available on our website.

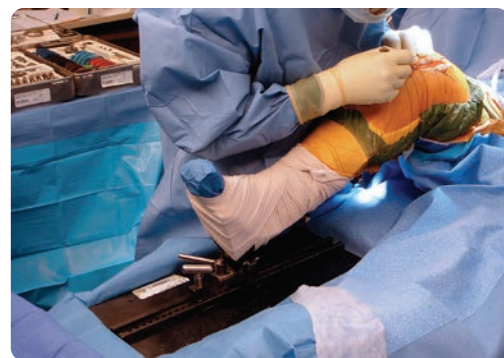


- ▶ Allows the leg to be manipulated into the desired position and securely locked in place
- ▶ Includes the necessary adjustments to tilt, rotate, and flex or extend the knee
- ▶ Extension/flexion adjustments can be made with the quick release of the ratchet
- ▶ In use, the base plate is clamped onto the operating table with the vertical side bar, the base plate is draped, and the sterile support plate lowered into the base plate
- ▶ The patient's foot is wrapped into the foot support with a sterile bandage (additional padding may be used for thin tibias)
- ▶ The complete unit is steam and gas sterilizable
- ▶ Three (3) Sterile Pads/Wraps are included with each new purchase (See below for Pad/Wrap information)

Tilt Bar

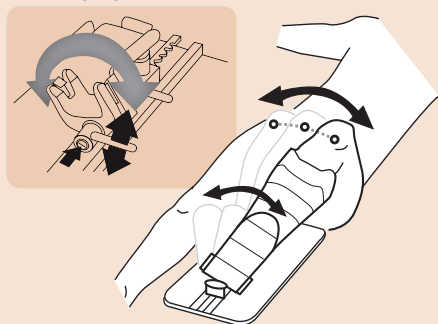


Loosening the Tilt Bar allows the knee to be tilted in either direction. Tightening the bar locks the Yoke System in the desired position.

Ratchet

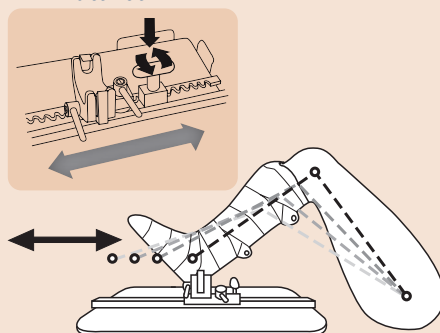

The Ratchet allows the Yoke Assembly to be moved in a precise gradual manner, the length of the Track. For faster adjustments, downward pressure on the Ratchet Handle releases the Yoke Assembly which then can be easily slid the length of the Track.

Rotation Bar

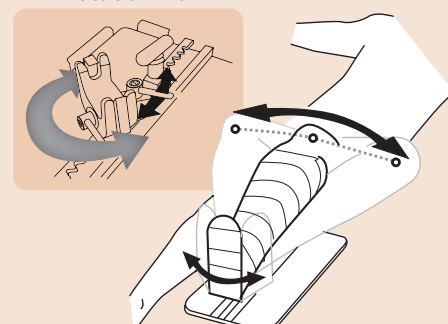

Loosening the Rotation Bar allows the knee to be rotated in either direction. Tightening the bar locks the Yoke System in the desired direction.

Leg Positioner Sterile Protective Pad & Wrap

Disposable, latex-free sterile foam pad and cohesive wrap helps protect patient from pressure sores, abrasions and possible neurological impairment while securing foot into the boot

PRODUCT NO:	
2629-00	[Case of 10 Sets]
2629-L	[1 Set]



Each set includes 1 Pad and 1 Wrap.



Compatible with the Innomed's Stulberg, Robb, Connelly, and TKR Leg Positioners





Nicholson Headrest Set

Designed by Gregory Nicholson, MD

Helps position the patient for all types of shoulder surgery in the beachchair position

Designed to provide excellent exposure to the shoulder, the headrest can be used with standard OR tables (with no modifications to the table). The headrest provides patient support and helps position the patient for all types of shoulder surgery—arthroscopic and open—in the beachchair position. It can be quickly placed and adjusted.



A gel pad forehead strap with velcro is included for optional use.

PRODUCT NO'S:

2450 [Headrest]

Main Plate Dimensions: 6" x 18" (15,2 x 45,7 cm)



Neck Offset Adjustment: 8" (20,3 cm)

Includes/Replacement Parts:

2450-S [Strap with Gel Pad]

4150-PD2 [Set of 2 Pads]

Set includes one Headrest, one Strap with Gel Pad, and two Pads.

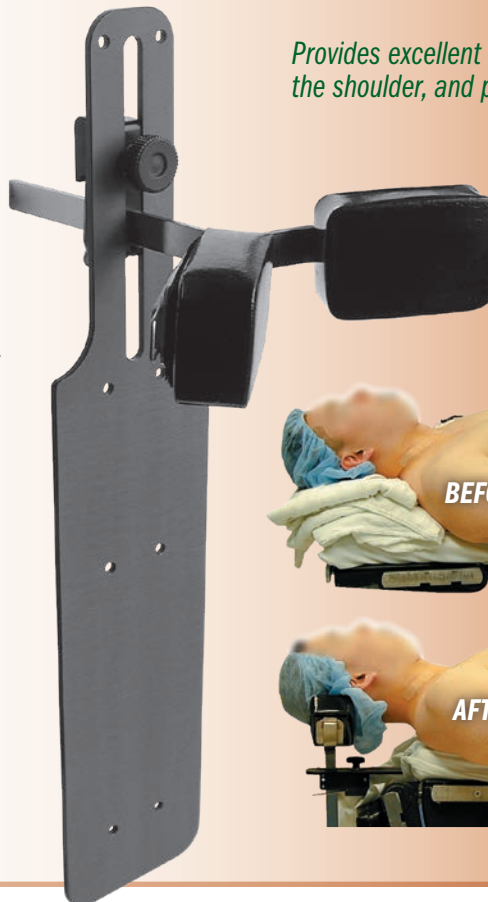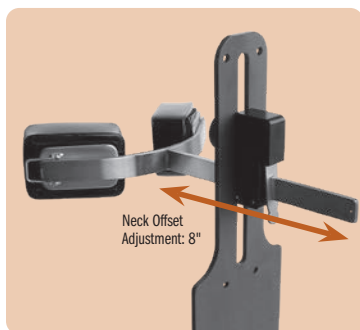

Provides excellent exposure to the shoulder, and patient support



Auerbach Arm Holder Rake Retractor Set

Designed by David M. Auerbach, MD

Allows intraoperative positioning for procedures of the posterior arm, elbow, and forearm



PRODUCT NO'S:

2415-00 [Set]

Individual/Replacement Parts:

2415-01 [Arm Holder Assembly]

Overall Length: 20" (50,4 cm)

Arm Holder Dimensions: 14.5" x 4" (36,9 x 10,2 cm)

Overall Width including Cleats: 7.5" (19,1 cm)

2415-02 [Arm Holder Upright Rod]

Overall Length: 19.25" (49,9 cm)

2415-04 [Rake Chain Retractor 4-Prong]

Two included in set, one with this product number

Overall Length including Chain: 10" (25,4 cm)

Retractor Width: .75" (1,9 cm)

2415-06 [Rake Chain Retractor 6-Prong]

Two included in set, one with this product number

Overall Length including Chain: 10" (25,4 cm)

Retractor Width: 1.25" (3,2 cm)


2595 [Table Clamp]

2770-P [Silicone Pad]

Dimensions: 12" x 5.5" (30,5 x 14 cm)

Replacement Parts:

2590-S [Black Straps] Pkg of 10

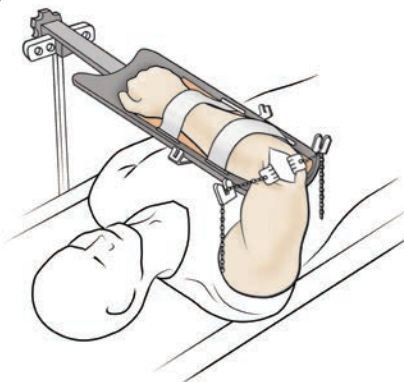

- ▶ Simple design for fast and easy positioning
- ▶ Connects over the drape in the sterile field using the supplied rail clamp and post
- ▶ Can be repositioned during surgery
- ▶ Sterilizable rubber pad protects the arm
- ▶ Retractors for the skin and soft tissues connect to the holder
- ▶ Compact for easy storage



Set includes: Arm Holder Upright and Assembly, (2) 4-Prong Chain Retractors, (2) 6-Prong Chain Retractors, (2) Black Straps, Silicone Pad, and Table Clamp.

Arm Board

A radiolucent arm board designed for all upper extremity surgery


- ▶ Securely slides under the table pad
- ▶ Patients body weight holds arm board in place

New!

PRODUCT NO:

2419

Overall Length: 42" (106,7 cm)
Maximum Width: 19" (48,25 cm)
Minimum Width: 10" (25,4 cm)
Thickness: .5" (1,27 cm)

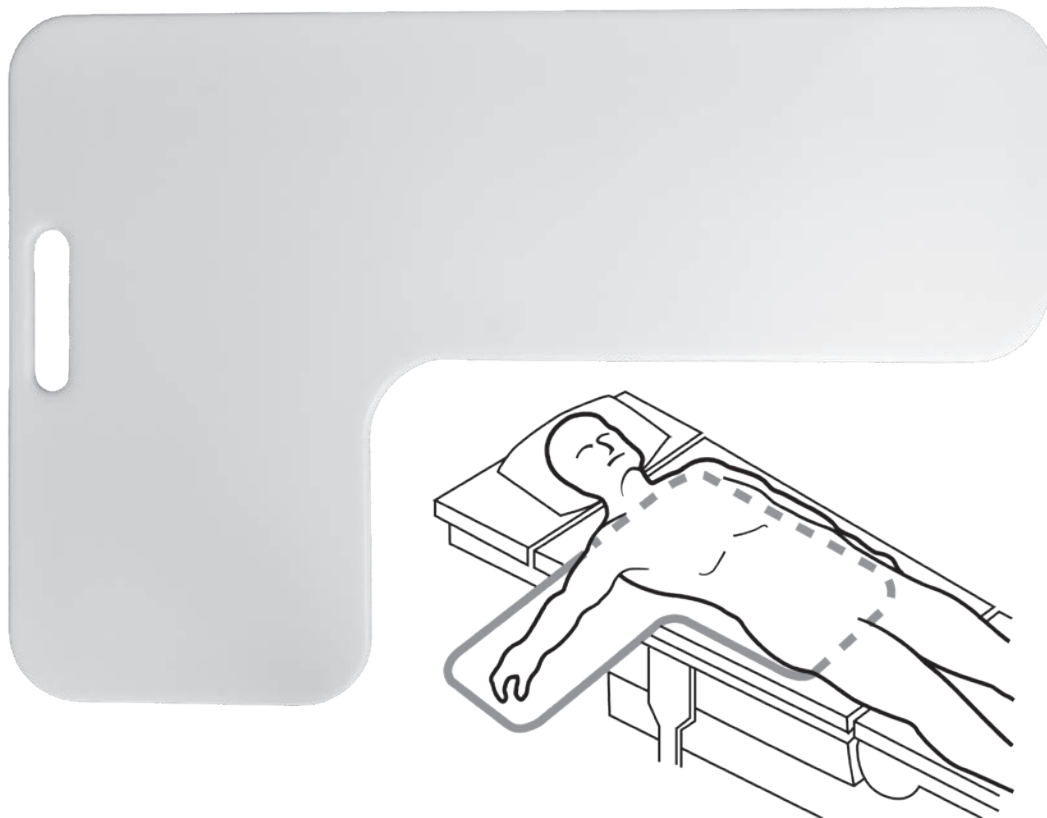

Distal Humerus Fracture Board

Designed by Burk Young, MD

Allows the surgeon to pin these fractures without having to manually hold the fracture reduced, allowing the surgeon to focus on accurate pin placement and reduction. The height of the crossbar is fully adjustable to accommodate different size patients. Reduction is achieved by an assistant gently applying axial traction through the forearm, with the crossbar applying the counter traction. Pinning is done with the C-arm in the lateral position. An optional separate attachment to support the arm for distal humerus fractures in adults is available. Unit not sterilizable.

PRODUCT NO'S:

2445 [Fracture Board – Pediatric]

Main Board Dimensions: 22" x 12" (55,8 cm x 30,5 cm)

Crossbar Height Adjusts From: 4.5" to 7.5" (11,4 cm x 19,1 cm)

2445-01 [Fracture Board – With Adult Adapter]

Optional/Replacement Part:



2445-06 [Adult Adapter]




Designed for the pinning of pediatric supra-condylar and adult distal humerus fractures



Board, crossbar, and optional adult adapter are radiolucent.

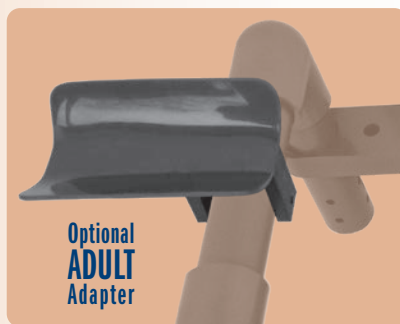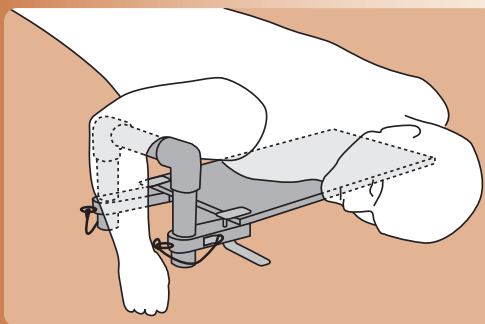

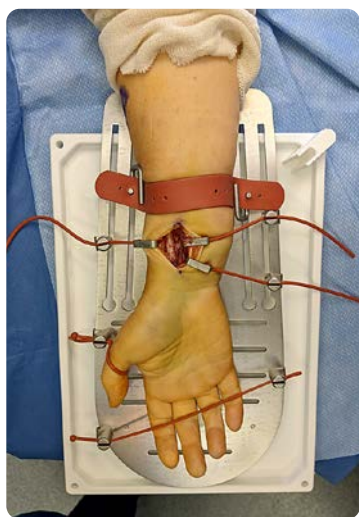

Thumb Post & Clip
Shown attached to plate


Suction Holder
Insert in any corner
to help remove blood
accumulating in tray



Auerbach Hand Positioner Set

Designed by David Auerbach MD

Designed to position and retract the skin for surgical exposures of the hand, wrist and forearm

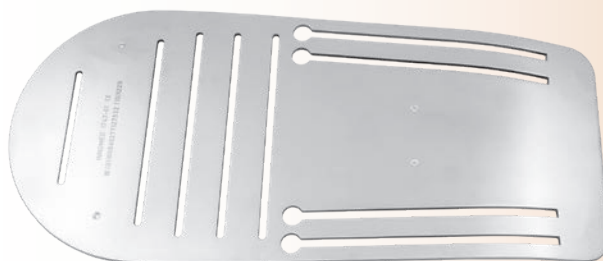

Hand Plate



Hand Tray



Thumb Post



Thumb Post Clip

Wrist Strap Buckles (2)



Suction Holder



Retractors (4)



Cord Clips (7)



Wrist Straps (2)



Cords (6)

PRODUCT NO'S:	
1747-00	[Auerbach Hand Positioner Set]
Also available individually:	
1747-01	[Hand Plate] Dimensions: 15" x 7" (38,1 x 17,8 cm)
1747-02	[Hand Tray] Dimensions: 13.75" x 9.75" (34,9 x 24,8 cm)
1747-03	[Thumb Post]
1747-03-C	[Thumb Post Clip]
1747-04	[Cord Clip] Seven (7) included in Set, One (1) with this product number.
1747-05	[Retractor] Four (4) included in Set, One (1) with this product number.
1747-06	[Wrist Strap Buckle] Two (2) included in Set, One (1) with this product number.
1747-07	[Wrist Strap] Two (2) included in Set, One (1) with this product number.
1747-08-6	[Set of 6 Cords]
1747-09	[Suction Holder]



Silicone Hand with Positioning Rings

Designed to help with positioning of hand and fingers for surgery, the silicone rings aid in stabilizing the fingers

The flexible silicone is easily bendable while maintaining the ability to remain in position once set. Silicone hand and rings are steam sterilizable.

PRODUCT NO'S:	
1746-00	[Silicone Hand with Pkg of 6 Rings]
Replacement Parts:	
1746-02	[Positioning Rings] Package of 12 Dimensions: 28 mm Outer Diameter, 3 mm Thick

MADE FOR INNOMED IN GERMANY



Quick Clip Rail Clamp

Designed for quick and easy attachment, the clamp snaps on to any OR side rail, and features a pull-to-release lever for quick adjustments



- ▶ Compatible with flat-blade style OR accessories and 1" x 1/4" blade post accessories
- ▶ Designed for patient weights up to 800 lbs.

PRODUCT NO:

2596

Dimensions:

3" Wide x 2" High x 2.75" Deep (max.)
7.6 cm Wide x 5.1 cm High x 7 cm Deep (max.)



New!

The GOAT Sterilizable Drape Clamp

Designed by William E. Nordt, III, MD

A sterilizable clamp that attaches to the OR table — designed to help eliminate the frustration and potential contamination of shifting drapes during surgery

PRODUCT NO:

2594-00 [Drape Clamp]

Dimensions with Screw:
2.5" Wide x 1.75" Tall x 2" Deep
(6.35 x 4.4 x 5.1 cm)

Includes / Available Individually:

2594-01 [Clamp Only]

Dimensions:
2.5" Wide x 1.75" Tall x 1.07" Deep
(6.35 x 4.4 x 2.7 cm)

4150-PS [Post Screw]

T-Handle Width: 1.875" (4.7 cm)
T-Handle Depth: 1.5" (3.8 cm)



New!

Table Clamps

Designed to help clamp and hold a device to the table

For Use with these Innomed Positioning Devices:

- ▶ Stulberg Leg Positioner
- ▶ Stulberg Hip Positioner
- ▶ Auerbach Arm Holder Rake Retractor Set
- ▶ Freeman Arm Holder
- ▶ Kirschenbaum Foot Positioners
- ▶ Robb Leg Positioner
- ▶ Thornberry Large Patient Hip Positioner

PRODUCT NO:

2595 [Table Clamp]



For Use with these Innomed Positioning Devices:

- ▶ Wixson Anterior Suspension Hook System
- ▶ Chandran Thigh Lift Positioner

PRODUCT NO:

9125 [Rotating Table Clamp]



For Use with these Innomed Positioning Devices:

- ▶ Capello Patient Positioner
- ▶ Direct Anterior Total Hip Arthroplasty Leg Positioner
- ▶ Durham Leg Positioner
- ▶ Leg Stabilizer
- ▶ Modified 90° Leg Stabilizer

PRODUCT NO:

9120 [Table Clamp]


9120-EU* [Table Clamp]



*9120-EU version fits European bed rails.

PRODUCT NO:

9121 [Clark Style Table Clamp]

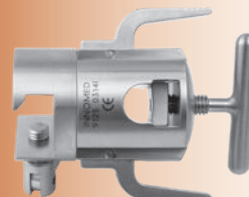

Reference Measurement Statement (Non-Measuring Function)

This device is not intended to provide quantitative measurements and does not incorporate a measuring function as defined under Regulation (EU) 2017/745 (MDR). Any scales, markings, or visual indicators present on the device are provided solely for reference and orientation purposes to support general use. These indicators are not calibrated against traceable standards and have not been validated for accuracy, precision, or repeatability.


Accordingly, values derived from these markings shall be considered approximate only and must not be used where accurate or traceable measurements are required for clinical decision-making, diagnosis, or other applications requiring defined metrological performance.

Users shall rely only on devices specifically intended and validated for measurement when quantitative accuracy is required.

Measurements in this Catalog

All effort has been made to ensure the accuracy of the measurements listed in this catalog, however, some small differences may exist between actual and listed measurements.

Measurements of **overall length** are the linear distance from one end of the product to the furthest opposite end, as shown in these examples:

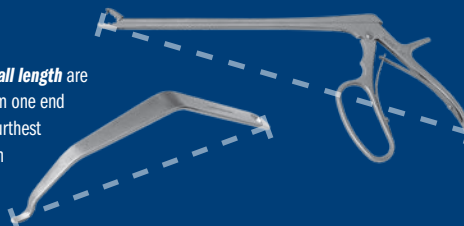

Measurements of **blade width** are the linear distance from one side of the product to the opposite side, typically at the widest point, as shown in this example:





103 Estus Drive
Savannah, GA 31404

PRSR.T. STD.
U.S. POSTAGE
PAID
BOURBON, IN
PERMIT NO. 29

FREE TRIAL

on most instruments

Instruments are available for a no-charge two-week evaluation – includes FREE UPS Ground Shipping*

*When shipped to a hospital or medical center; additional charge applies for expedited shipping.
Free trial offer excludes implant extraction instruments, which are available as rentals. There is a pad replacement charge with the hip positioners.

New!

TKR Knee Positioner

Provides stable positioning of the knee during surgery

Three (3) Sterile Pads/Wraps are included with each new purchase.
(See page 15 for Pad/Wrap information)



PRODUCT NO'S:	
2618-00	[Set]
Set Includes:	
2618-FP	[Aluminum Pivot Mount Footpiece] Length: 12.93" (32,8 cm) Height: 8.61" (21,9 cm) Width: 6.2" (15,7 cm)
2618-BP	[TKR Leg Positioner Base Plate] Length: 24" (61 cm) Height: .875" (2,22 cm) Width: 6.55" (16,6 cm)
2618-01	[TKR Foot Clamp Slide] Length: 3.75" (9,5 cm) Height: 4" (10,2 cm) Width: 4.375" (11,1 cm)
2629-L	[Sterile Pad/Wrap] Three included in set; product number is for one only
Optional & Replacement Parts:	
2629-00	[Case of 10 Sterile Pads/Wraps]



ISO 13485:2016

INNOMED, INC
103 Estus Drive
Savannah, GA 31404

Tel 912.236.0000
Fax 912.236.7766

innomed.net
info@innomed.net

Innomed-Europe LLC
Alte Steinhauserstr. 19
CH-6330 Cham, Switzerland
Tel 0041 (0) 41 740 67 74
www.innomed-europe.com
orders@innomed-europe.com

Innomed-Europe GmbH
c/o Emons Logistik GmbH
In Rammelswiesen 9
D-78056 Villingen-Schwenningen, Germany
Tel 0049 (0) 7720 46110 60
www.innomed-europe.com
orders@innomed-europe.com

07
26

© 2026 Innomed, Inc.
All Rights Reserved

1.800.548.2362



INNOMED.NET